

EXHIBIT EE

ETHICON, INC.

a Johnson & Johnson company

P.O. BOX 151
SOMERVILLE • NEW JERSEY • 08876-0151

October 15, 1992

cc: B. Matlaga
J. McDivitt
↓
A. Melveger
RDCF

Mark Cafone

SEVEN YEAR DATA FOR TEN YEAR PROLENE™ STUDY: ERF 85-219

This report contains a summary of IR, IV, GPC, OM and SEM data supporting this study.

IR and IR Microspectroscopy (D.Burkley)

IR examinations were done for all explants at all sites to verify the suture identity for each explant. For all explanted sutures recovered from all 6 sites for every dog in this study, IR data showed each suture to be correctly identified.

IR microspectroscopy was used to examine cracked areas in ETHILON, Novafil and PROLENE™ explants. IR spectra obtained for cracked PROLENE specimens (Figure A) showed possible evidence of slight oxidation (a broadened weak absorbance at about 1650 cm⁻¹). IR spectra obtained for cracked areas of ETHILON and Novafil did not differ from uncracked areas (Figures B and C), but expected IR absorbances for oxidation would be masked by the strong carbonyl absorbances normally observed for these sutures. Figures D and E show pictures of the areas examined by IR microspectroscopy for ETHILON and Novafil.

IV and GPC (E.Muse)

Gel Permeation Chromatography (GPC) was run on PROLENE sutures explanted from dogs after seven years. The GPC data was compared to data from a current 4/0 PROLENE suture. The results indicate that there was no significant difference in molecular weight between the 4/0 PROLENE control and the seven year explants.

The following PROLENE explant samples were analyzed:

- Dog 1995 - site 3 (SR33853)
- Dog 2007 - sites 1 and 6 (SR34003)
- Dog 2008 - site 2 (SR34066)
- Dog 2019 - sites 2 and 3 (SR34180)

The GPC analysis was run on the Waters 150C GPC at 140°C using 1,2,4 trichlorobenzene as a mobile phase with Waters GPC columns. The instrument was calibrated with polypropylene standards.

Inherent Viscosity (IV) was determined on ETHILON™ and Novafil sutures explanted from dogs after seven years. The IV data¹ was compared to IV data from one and two year explants. The following results were found:

- 1) No significant differences were seen in IV values after one and two years.
- 2) Seven year IV values ranged from 75% to 93% of the one and two year IV values for ETHILON sutures.
- 3) Seven year IV values ranged from 75% to 90% of the one and two year values for Novafil.

The dog explant samples examined were from duplicate sites on four dogs for each time period (one, two and seven years). The IV data was determined using concentrations of 0.1 dl/g with HFIP as a solvent at 25°C.

OPTICAL MICROSCOPY and SCANNING ELECTRON MICROSCOPY (E.Lindemann)

Conclusions

- The 7 year in-vivo results generally substantiated the five year findings. They also closely correspond to the observations of explanted sutures from the dog that died prematurely after 6 years and 10.5 month implantation time.
- Degradation in PROLENE is still increasing and PVDF, even though a few cracks were found, is still by far the most surface resistant in-house made suture in terms of cracking.
- Of the eight explanted ETHILON sutures all showed heavy cracking and, in many cases, abrasion of the dyed surface layer. A decrease in the suture diameter was apparent in several cases.
- Cracks were not found in the seven Novafil explants. However a few longitudinal scratches probably due to mechanical damage and one longitudinal crack were observed.

Introduction

In November 1985 twenty-four dogs had been implanted with sets of ETHILON, PROLENE, PVDF and Novafil sutures for a ten year study. In 1990, after five years, explants from 5 beagle dogs were described in "TEN YEAR IN-VIVO STUDY SCANNING ELECTRON MICROSCOPY FIVE YEAR REPORT" by Elke Lindemann. The next explantation, after 7 years, was to start in June 1992. However, after 6 years and 10.5 months dog #1995 died prematurely. The microscopical examination of those explants was described in "TEN YEAR IN-VIVO STUDY: SCANNING ELECTRON AND LIGHT MICROSCOPY INTERIM REPORT ON DOG #1995 AFTER 6 YEARS, 10.5 MONTH, SR# 33788 and are included in the conclusion section of this report. In June of 1992 after 7 years, sutures were explanted from another set of 4 dogs. This report presents the results of the light and scanning electron microscopical examination of those explants.

¹SR33853, SR34003, SR34066, SR34180

Experimental

Four dogs had been implanted in November 1985 with the following 5-0 sutures:

	Site 1	Site 2	Site 3	Site 4	Site 5	Site 6
Dog 2001	PVDF	ETHILON	Novafil	PROLENE	PROLENE	Novafil
Dog 2007	PROLENE	Novafil	ETHILON	PVDF	PVDF	PROLENE
Dog 2019	Novafil	PROLENE	PROLENE	PVDF	ETHILON	ETHILON
Dog 2008	ETHILON	PROLENE	Novafil	PVDF	ETHILON	PVDF

Starting in June of this year the above dogs were sacrificed in weekly intervals. Approximately 20cm long sections of the explanted sutures were received in microscopy in glass vials which were kept refrigerated until they were examined.

Also the explanted LC 100 clip with about 2cm of each suture bundle was delivered in the same vial. The clip and the attached sutures were still deeply embedded in the surrounding tissue. These 'not cleaned' sutures were supposed to answer the question whether the process of cleaning and tissue removal might be responsible for an observed cracking. The primary concern of this study was however to examine the long pieces of explanted suture. Most of these specimens were still surrounded with some tissue, fortunately at a level low enough not to obscure examination in the light microscope under transmitted light. It was possible to examine the embedded PROLENE suture where the cracking of the suture was seen through the tissue. For this reason and time constrains the clip-attached sutures were not examined at this time.

To show that the drying and coating with a metal under vacuum, necessary for SEM examination, did not introduce cracking and other surface defects each strand of each long suture was 100% inspected in the Olympus Light Microscope in water. Oil, the usual medium for light microscopical inspection, was not chosen for this examination in order to eliminate surface changes during sample preparation. To cut down on lensing effects of the curved suture, the samples were photographed in polarized light using a 10x phase condenser with an ordinary transmitted light 20x objective (a 20x phase condenser was not available). The light diffraction introduced by the phase condenser was enough to allow an easier focusing at the focal plane of the largest diameter. Photomicrographs were prepared at 285x of areas which showed surface changes.

Strands of the suture including the above areas were then prepared for SEM observation in the JEOL JSM 840 AII by coating them under vacuum with gold to provide an electron conductive surface. Photomicrographs were prepared at 500x magnification.

Results

1) LM and SEM of PROLENE suture explants from seven implantation site.

In Figure 1A through 1D one area per site from each of the four dogs is shown in transmitted light. Out of seven sites cracking was found on PROLENE sutures from three sites. Notice the cracks observable through the still adhering tissue in Figure 1A in the suture from site 2.

In Figure 1 and 2 SEM views of areas are shown after most of the tissue had been carefully removed. Again out of seven sites sutures from three sites had areas which showed cracking.

2) LM and SEM of ETHILON suture explants from six implantation sites.

In Figure 3A through 3C sutures are shown from six different sites. Transmitted light allowed visualization of the differences between the intact dyed surface layer and the underlying colorless layers of the suture. In Figure 3A site 5 and Figure 3C site 3 the colorless area had not only lost its dyed surface layer but was abraded to such a degree that a decrease in suture diameter was found.

In Figures 3 and 4 the cracking and abrasion on sutures from all six sites, as observed with the SEM, is shown. Here also the decrease in diameter is particularly dramatic in Figure 3 site 1.

3) LM and SEM of PVDF suture explants from six implantation sites.

Figure 5A through 5C show six sites of PVDF explants as seen with the light microscope. Notice the intact surface on all the sutures.

In Figures 5 and 6 the SEM examination of the PVDF sutures is shown. Only on the suture from one site (Figure 6 site 6) some cracks are found. The surfaces of the sutures from the other five sites show some striations which could be mechanical damage, otherwise the surfaces look intact. The contaminant on the site 4 (Figure 5) suture is tissue which had not been removed completely.

4) LM and SEM of Novafil suture explants from five implantation sites.

Figure 7A through 7C show the Novafil sutures as observed with the light microscope. All surfaces from all sites look undamaged. Figure 7 and 8 show the SEM examination of these sutures. A few longitudinal scratches and cracks were found, see sites 1,2,3 (Figure 7,8). Also on the site 2 suture (Figure 8) still adhering tissue is found.

5) Degradation dependency on implantation site

To probe the question as to whether one implantation site might be more or less stressful towards the suture, a comparison was made of the six sites.

	Site 1	Site 2	Site 3	Site 4	Site 5	Site 6
Dog 1995	ETHILON cracks	PVDF	PROLENE cracks	Novafil	Novafil cracks	ETHILON cracks
Dog 2001	PVDF	ETHILON cracks	Novafil	PROLENE	PROLENE cracks	Novafil
Dog 2007	PROLENE	Novafil scratch	ETHILON cracks	PVDF	PVDF	PROLENE cracks
Dog 2019	Novafil scratch	PROLENE	PROLENE	PVDF	ETHILON cracks	ETHILON cracks
Dog 2008	ETHILON cracks	PROLENE cracks	Novafil cracks	PVDF	ETHILON cracks	PVDF cracks

The only site, in the 5 dogs of this study, from which sutures were explanted that showed no surface damage was site 4. However, of those five sutures three were PVDF and one was Novafil. Those are the sutures that showed only marginal surface changes in this study. Therefore this observation can be discounted.

Elke Lindemann

Eugene P. Muse

Daniel F. Burkley

Attachment

7 YEAR.DFB

Figure A - Comparison IR Spectra of Cracked and Non-cracked Regions of PROLENE (SR33788)

C:\FIRST\SH\WORK2\33788-2

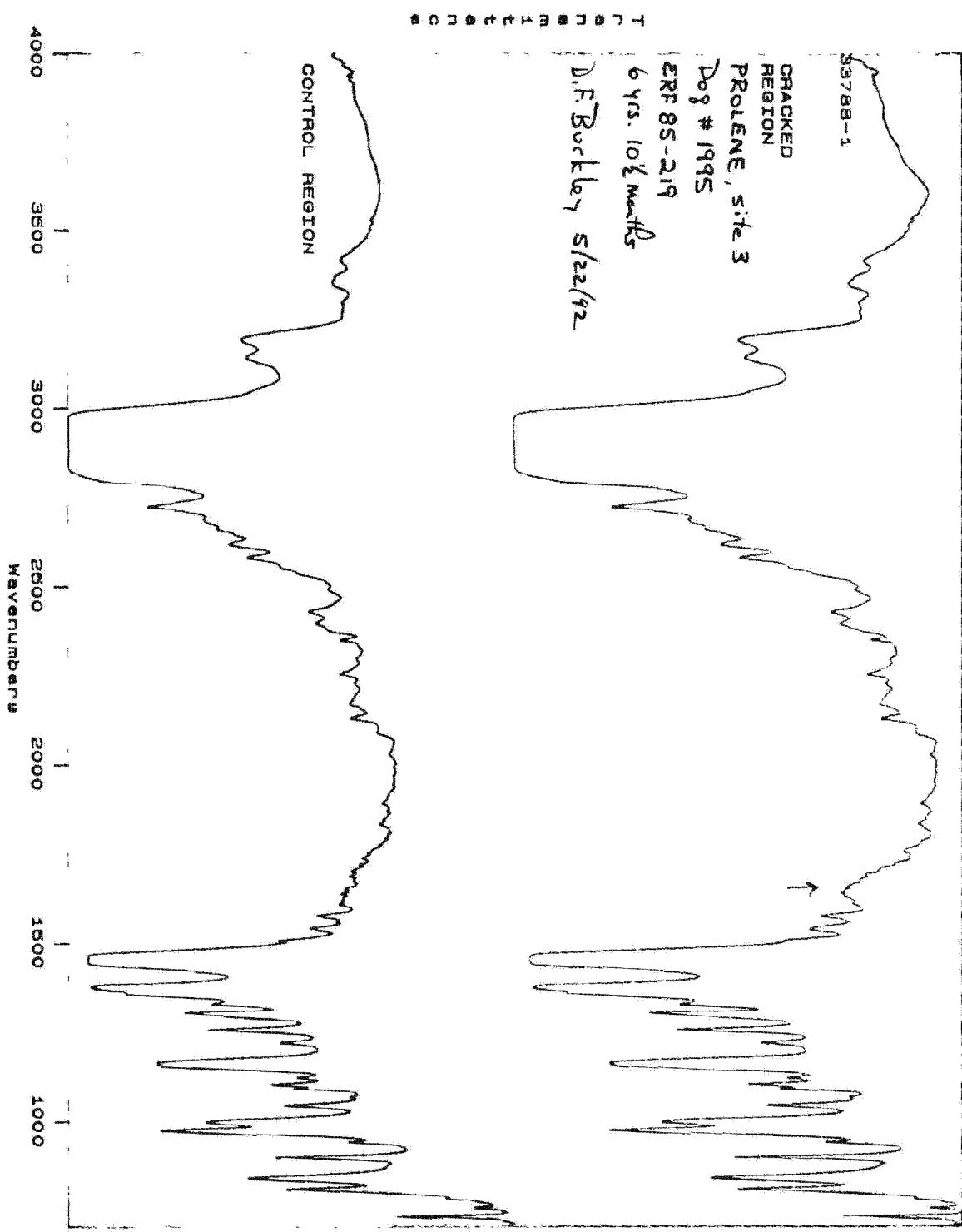


Figure B - Comparison IR Spectra of Cracked and Non-cracked Regions of ETHILON (SR33853)

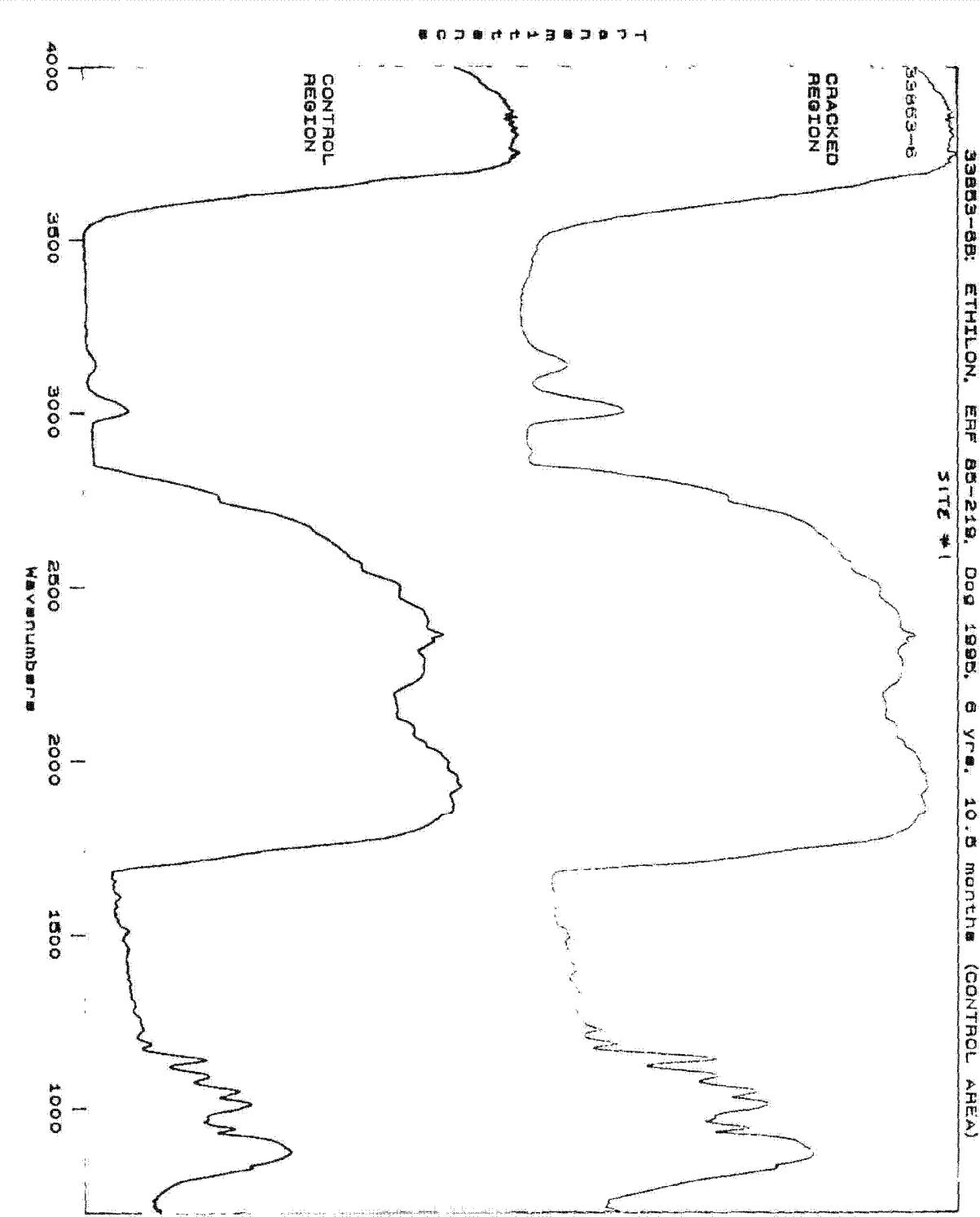
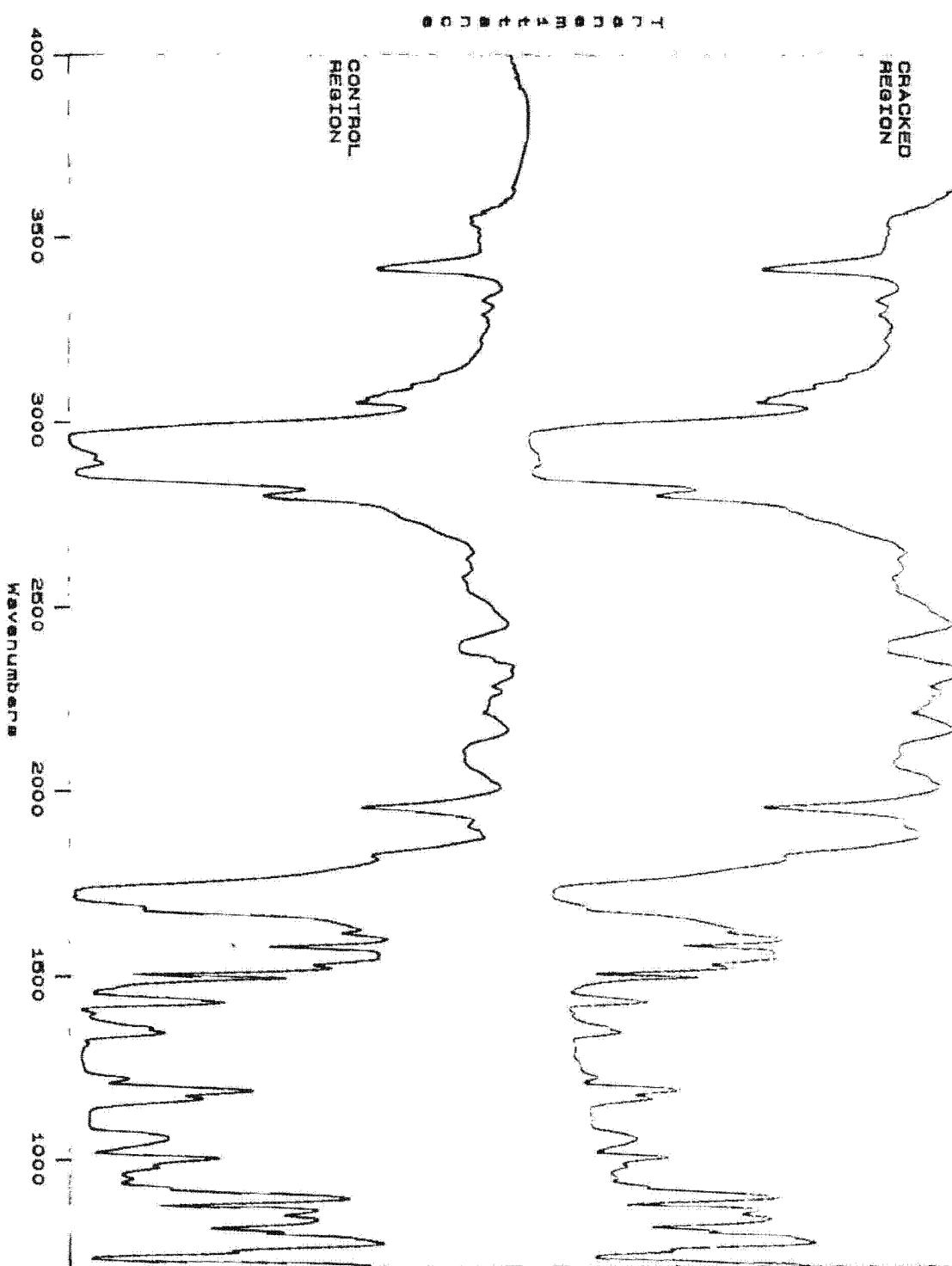


Figure C - Comparison IR Spectra of Cracked and Non-cracked Regions of Novafil (SR33853)

33853-1b: NOVAFIL. ERF 88-218. Dog 1995. 8 yrs., 10.5 months (CONTROL AREA)



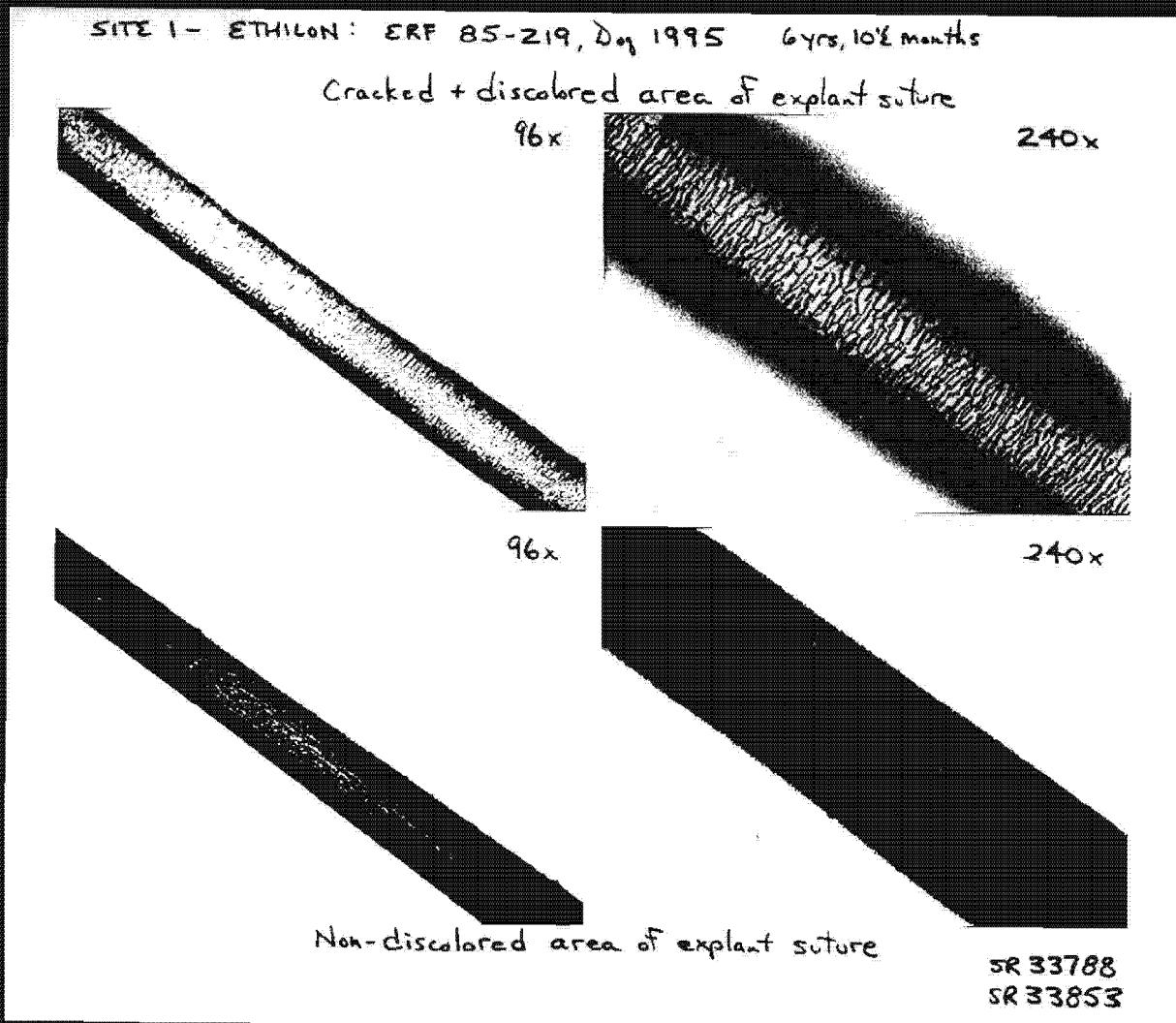


Figure D - Comparison Pictures of Cracked and Non-cracked Regions of ETHILON (SR33853)



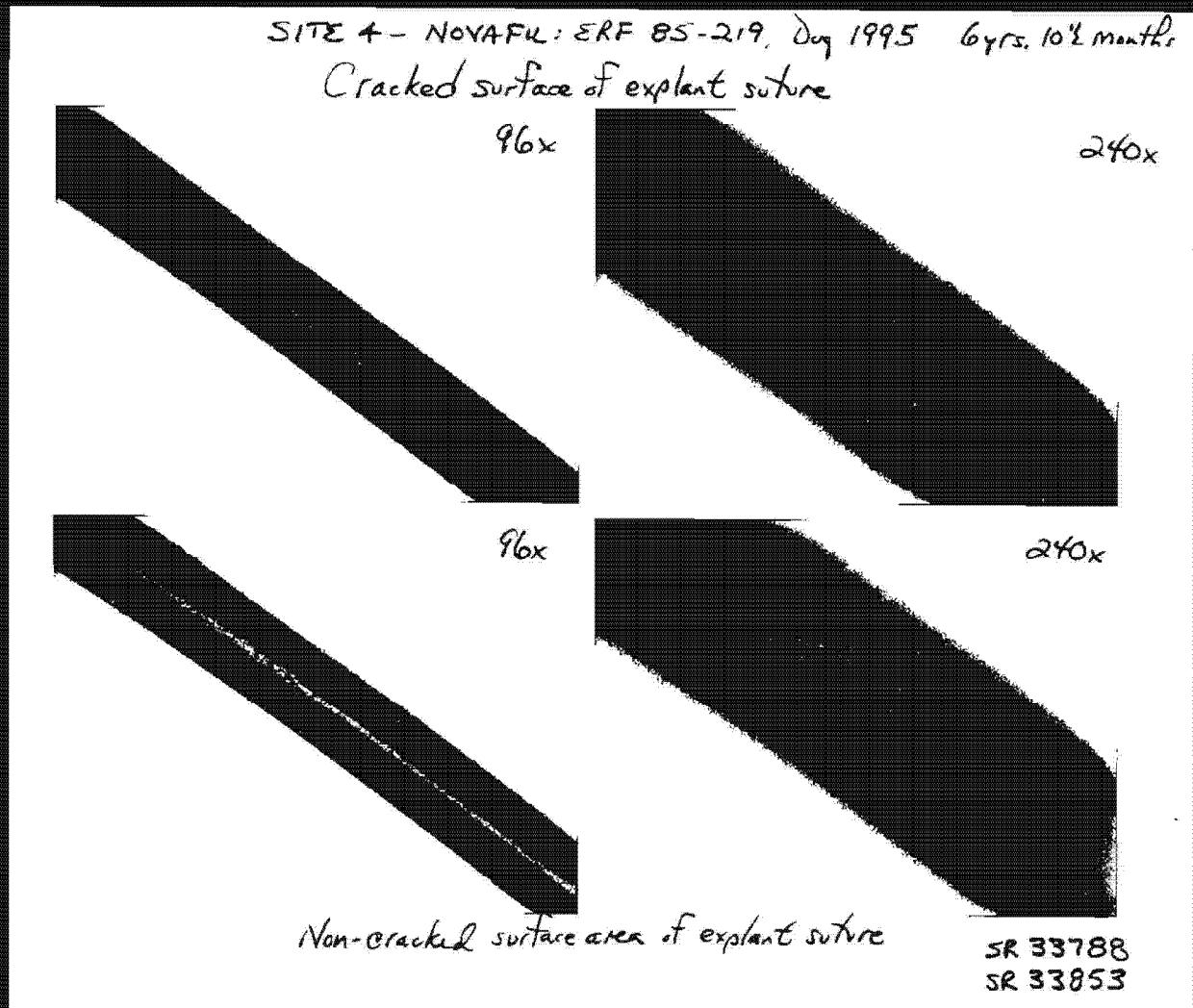


Figure E - Comparison Pictures of Cracked and Non-cracked Regions of Novafil (SR 33853)

Figure 1A 7 Year Prolene Explants

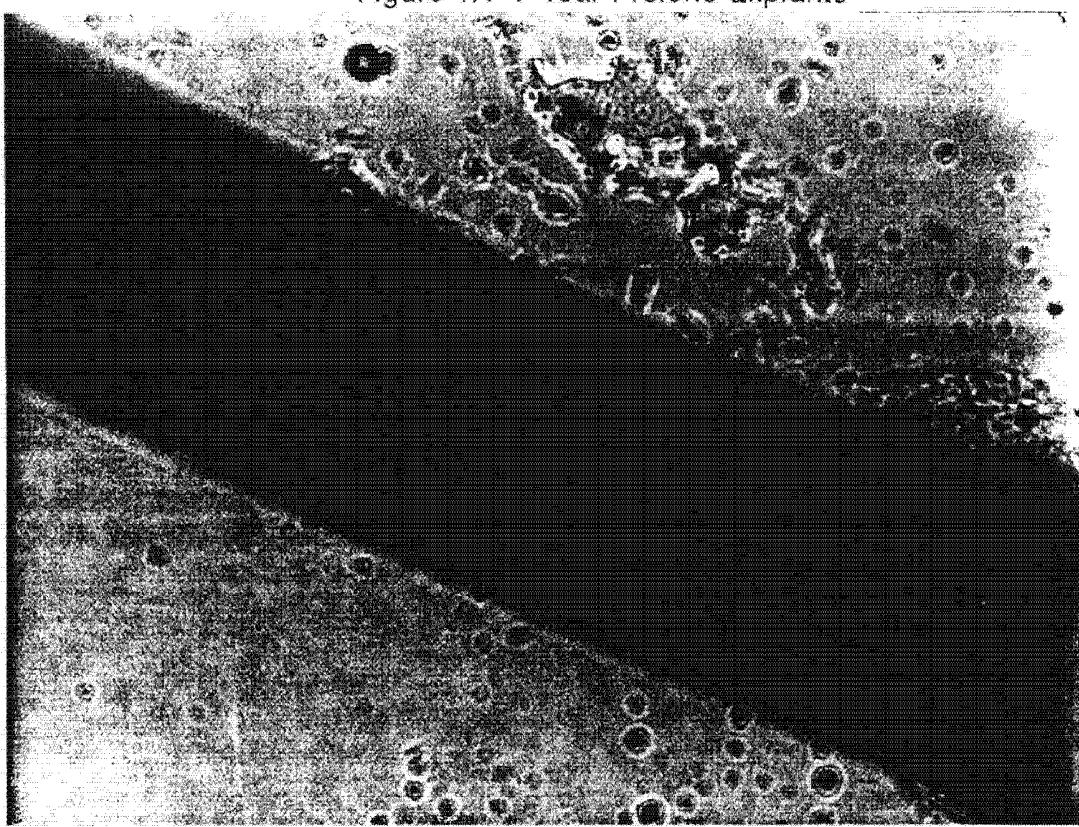

dog 2001 site 5

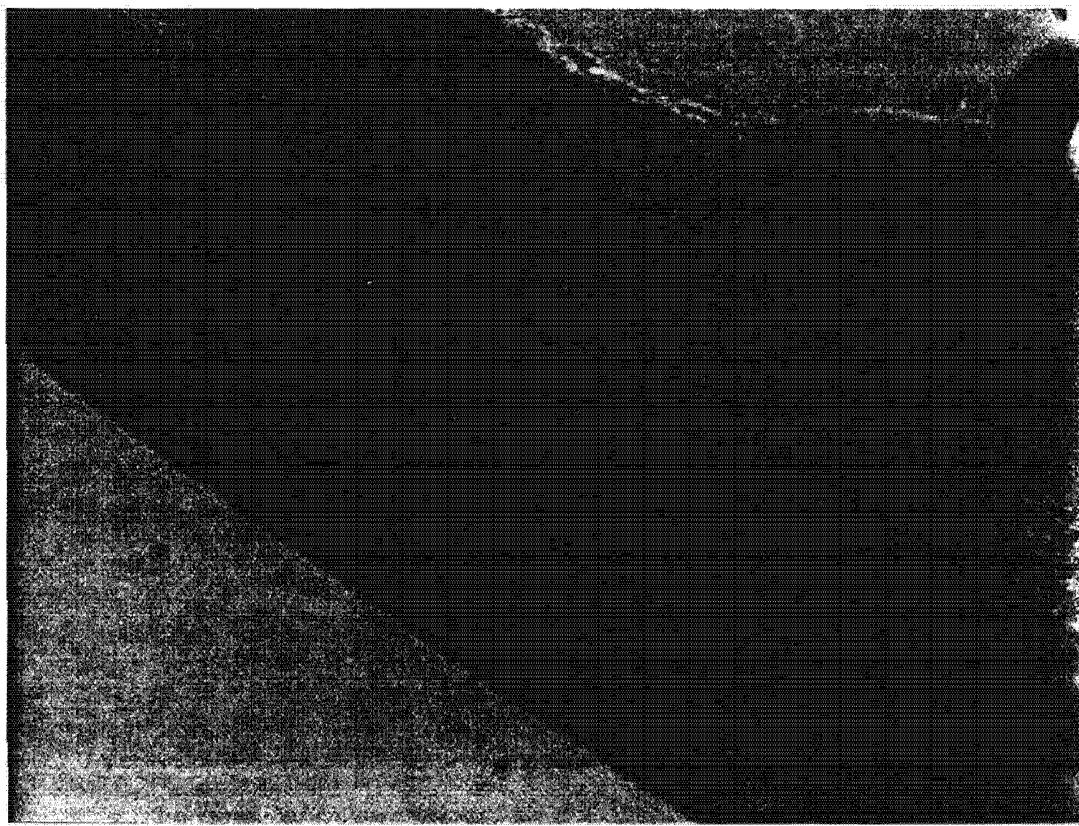


dog 2008 site 2

E. Lindemann 7/9/92 SR# 33985

CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690762

Figure 1B 7 Year Prolene Explants

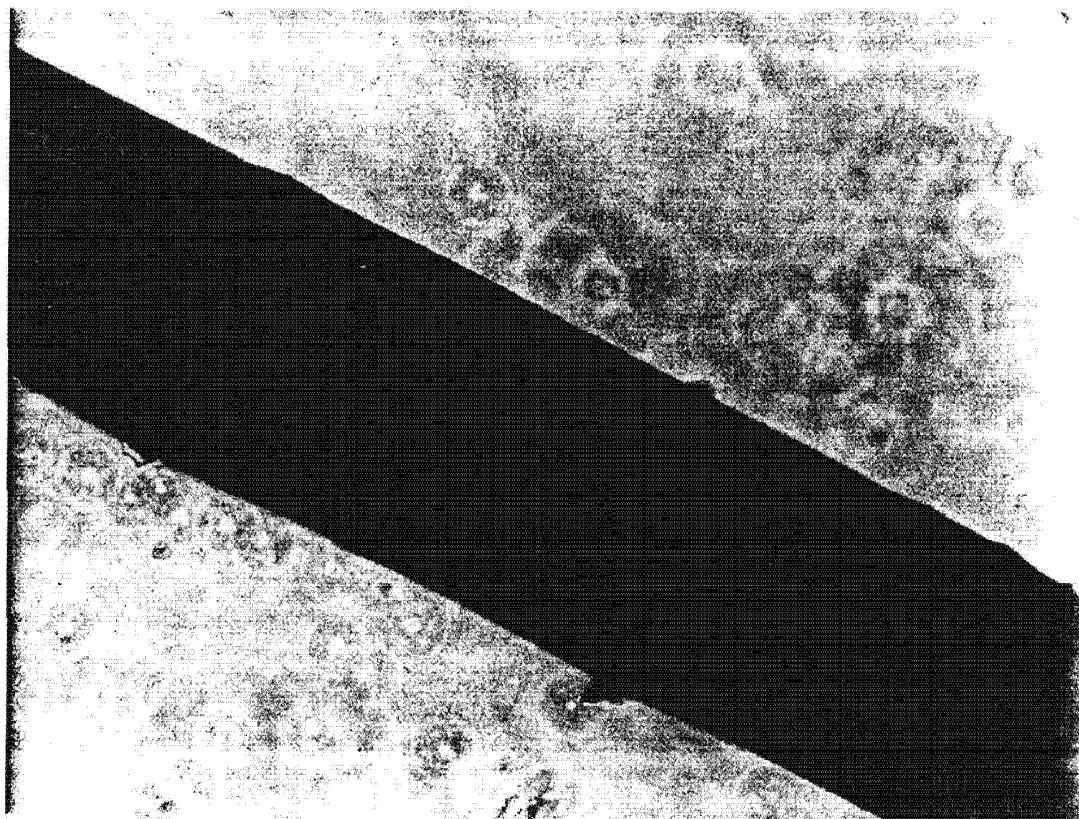


E. Lindemann 7/9/92 SR# 33985


CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690763

Figure 1C 7 Year Prolene Explants

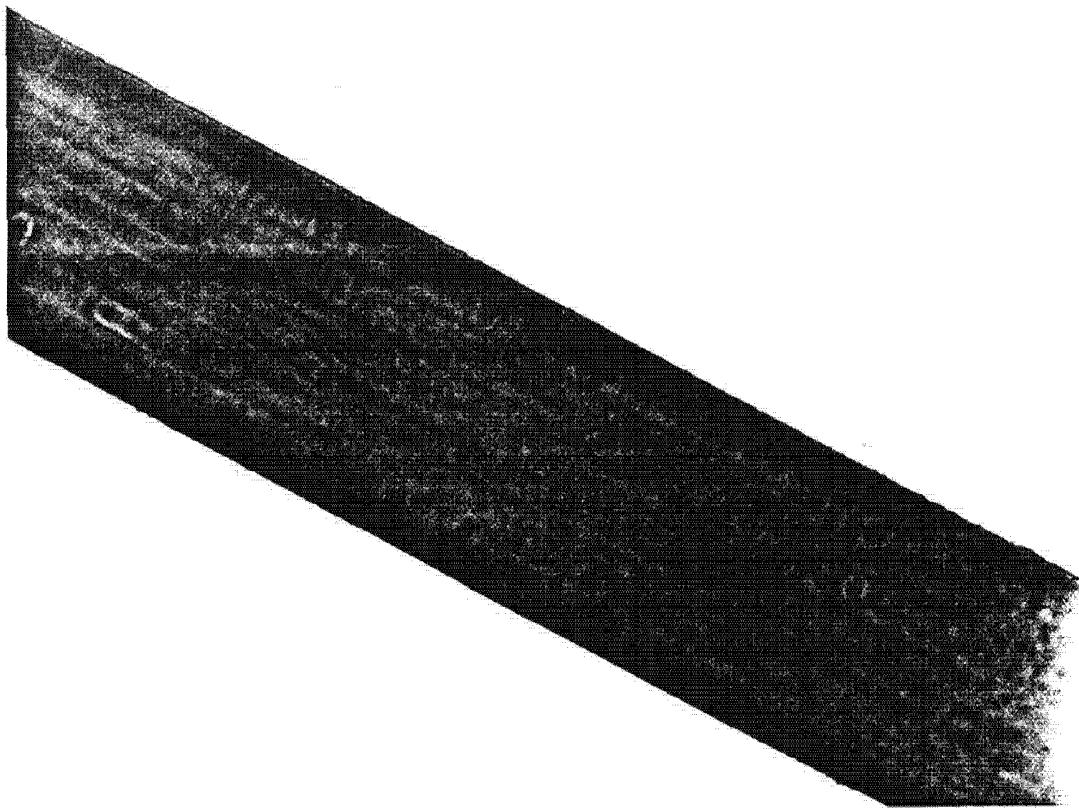

dog 2007 site 1

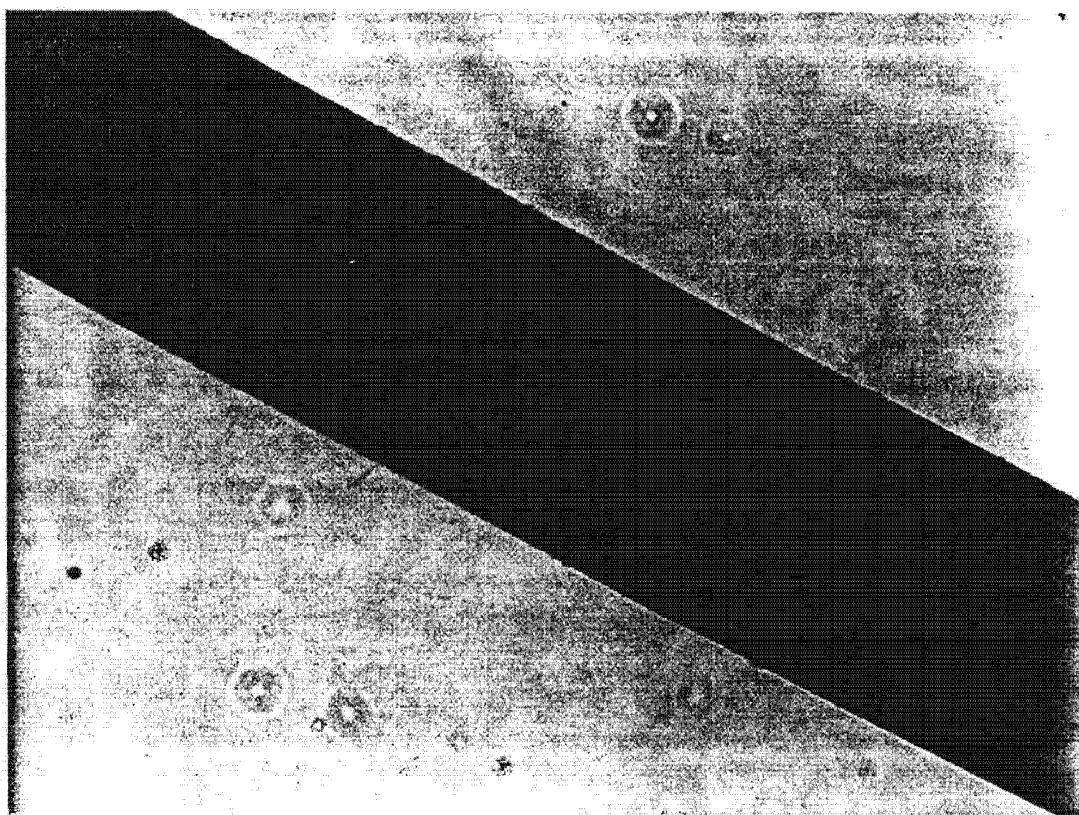

dog 2019 site 3


E. Lindemann 7/9/92 SR# 33985

CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690764


Figure 1D 7 Year Prolene Explants



dog 2001 site 4

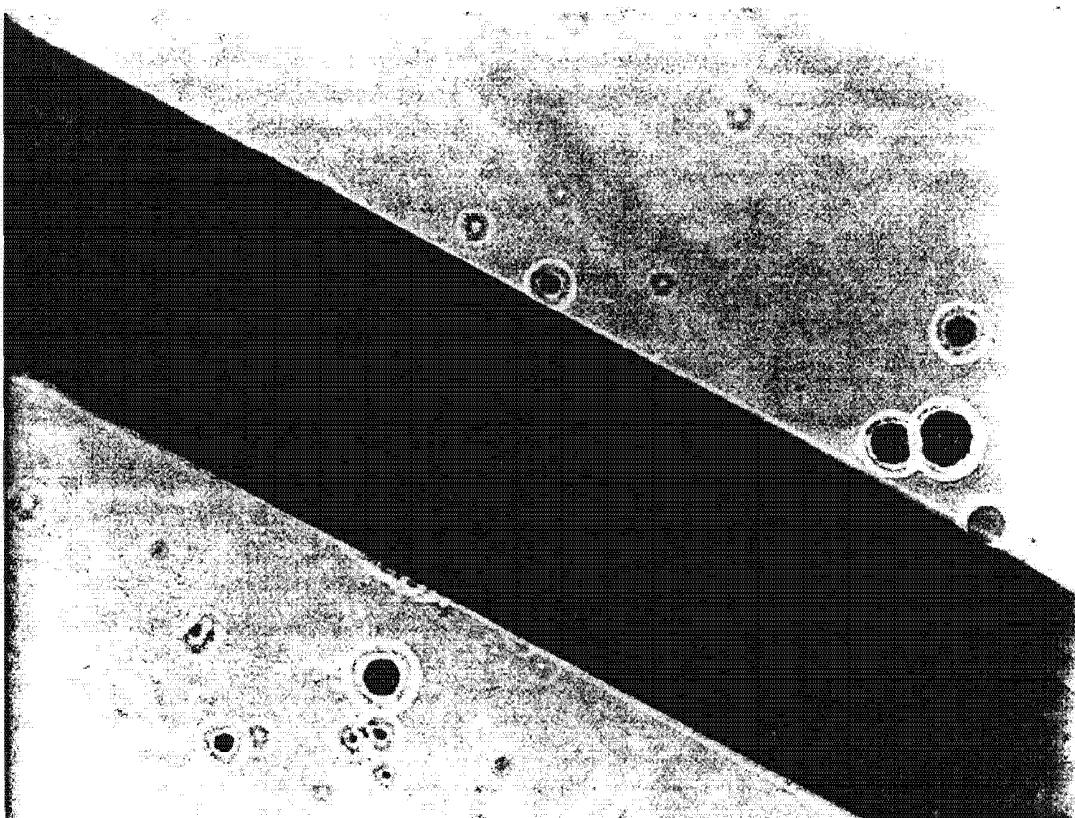



E. Lindemann 7/9/92 SR# 33985

Figure 1 7 Year Prolene Explants

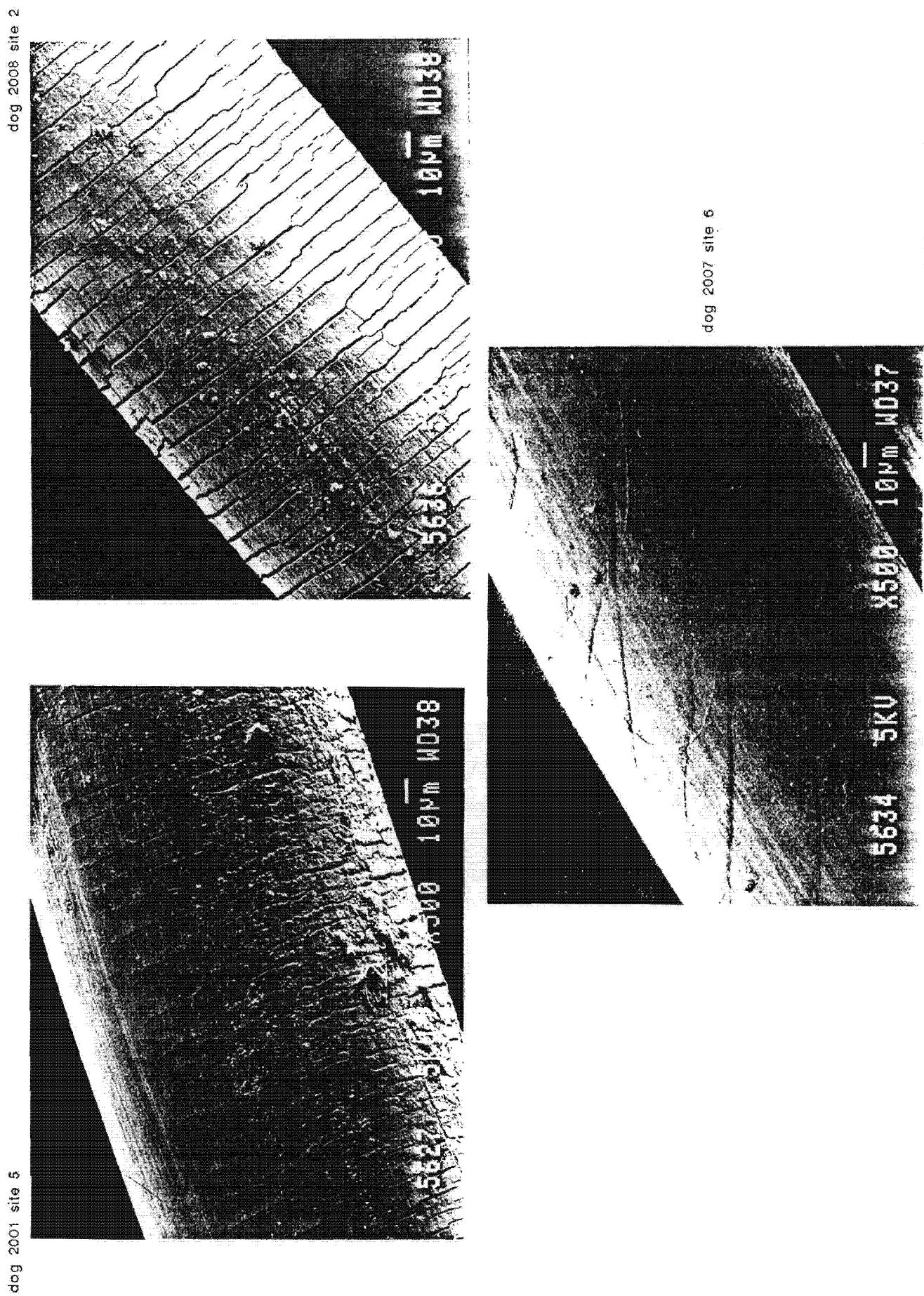

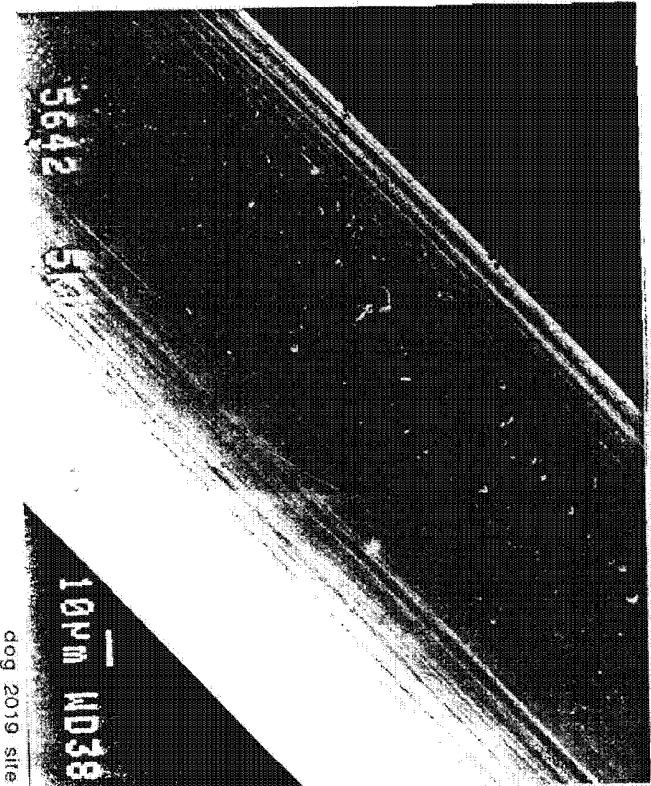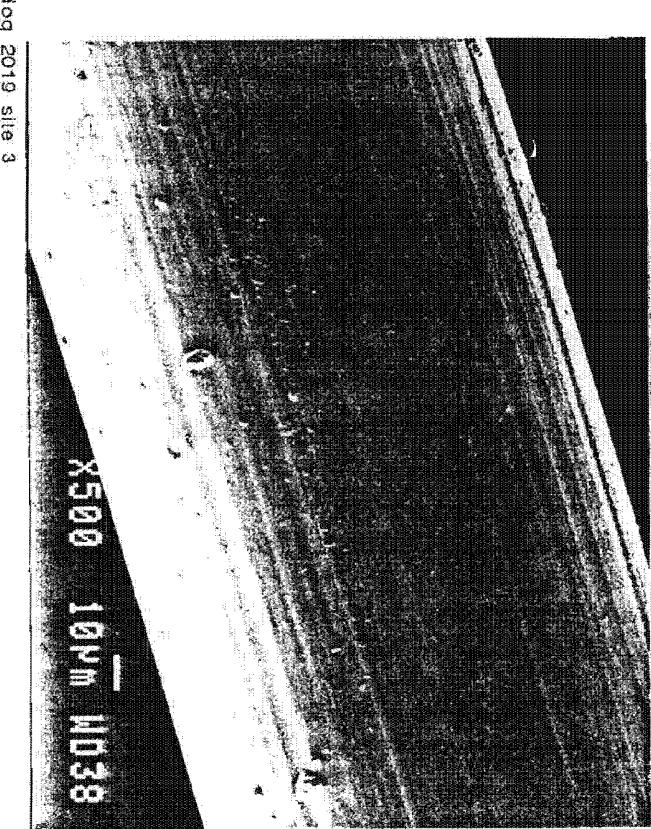


dog 2019 site 2

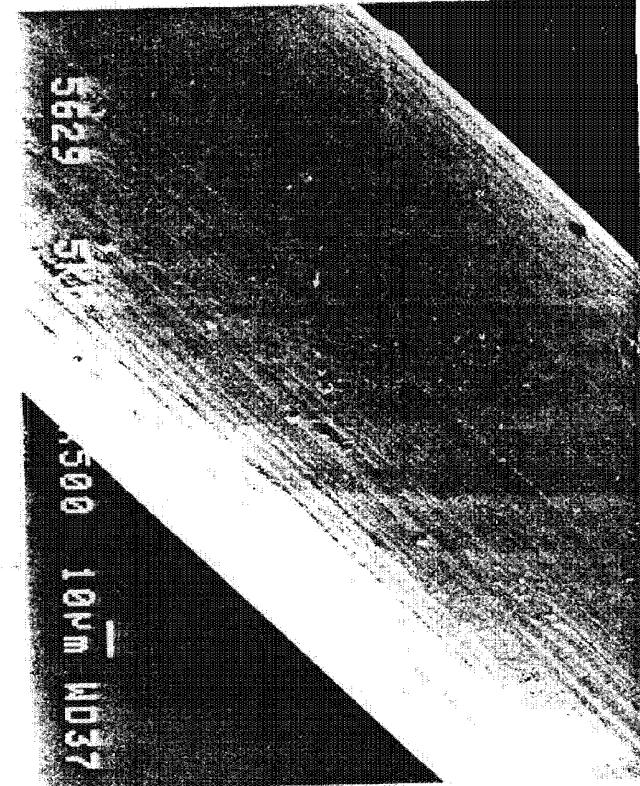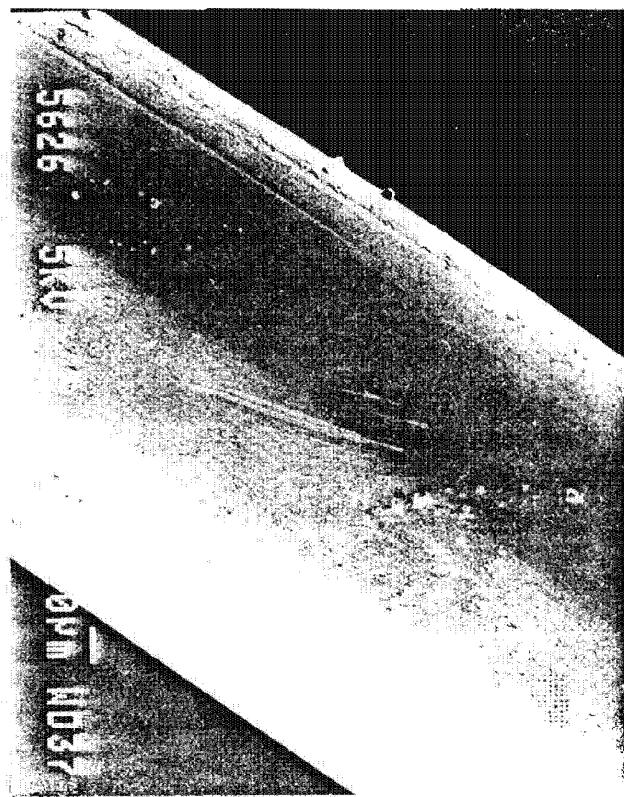


E. Lindemann 7/9/92 SR# 33985
dog 2007 site 1

Figure 2 7 Year Prolene Explants

Figure 3A 7 Year Ethilon Explants

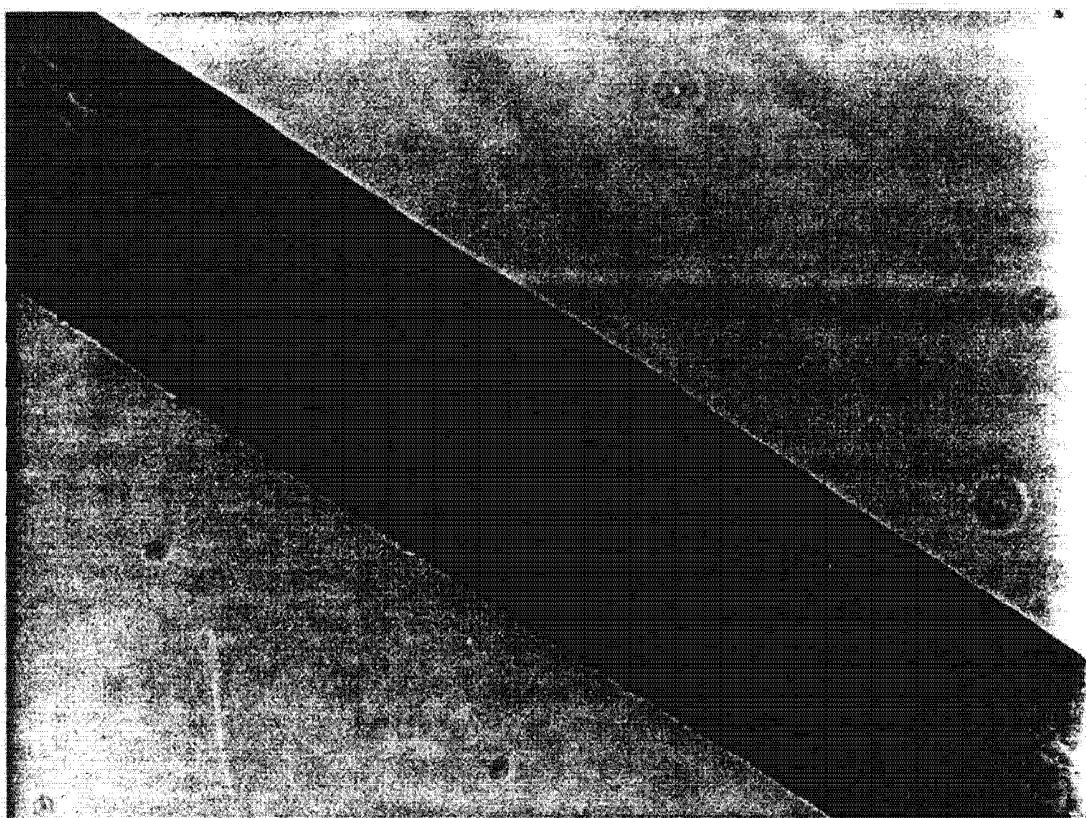

dog 2019 site 6

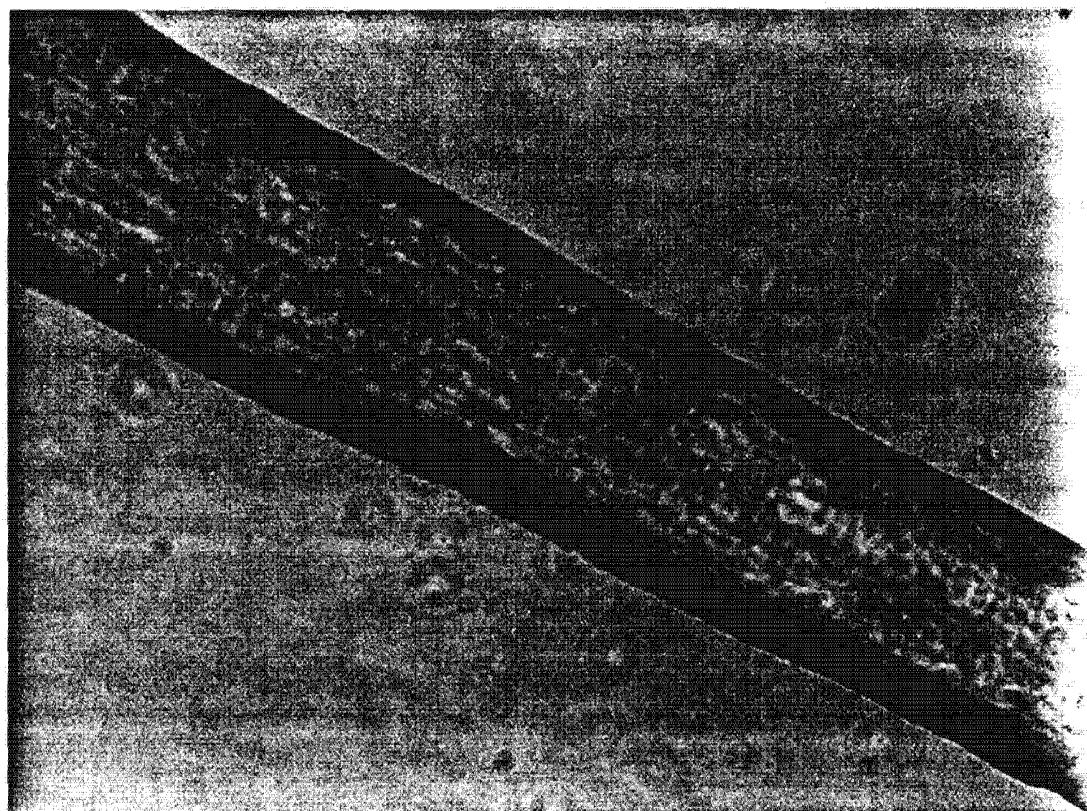


dog 2008 site 5

E.Lindemann 7/15/92 SR# 34097

CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690768

Figure 3B 7 Year Ethilon Explants

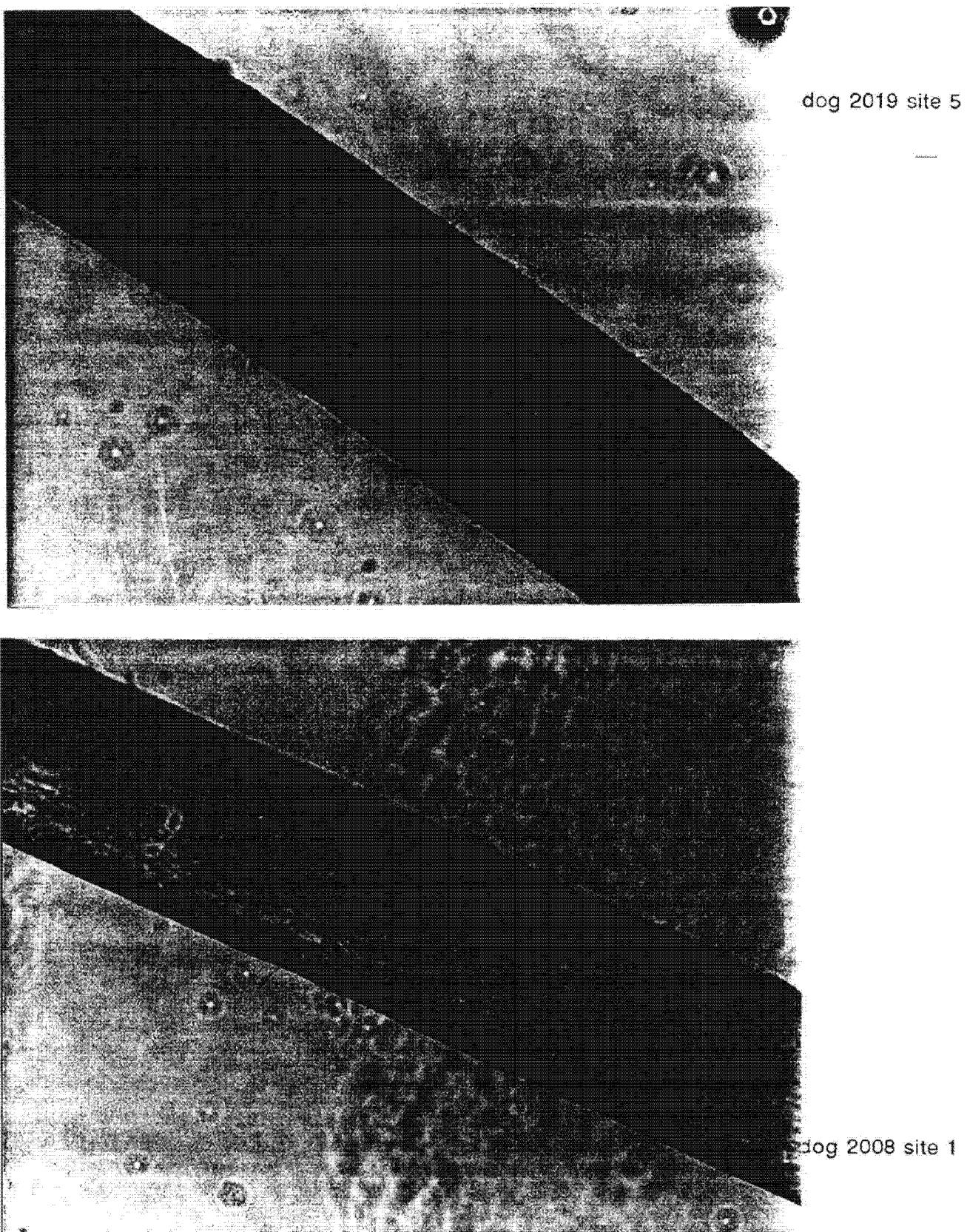


E.Lindemann 7/15/92 SR# 34097


CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690769

Figure 3C 7 Year Ethilon Explants

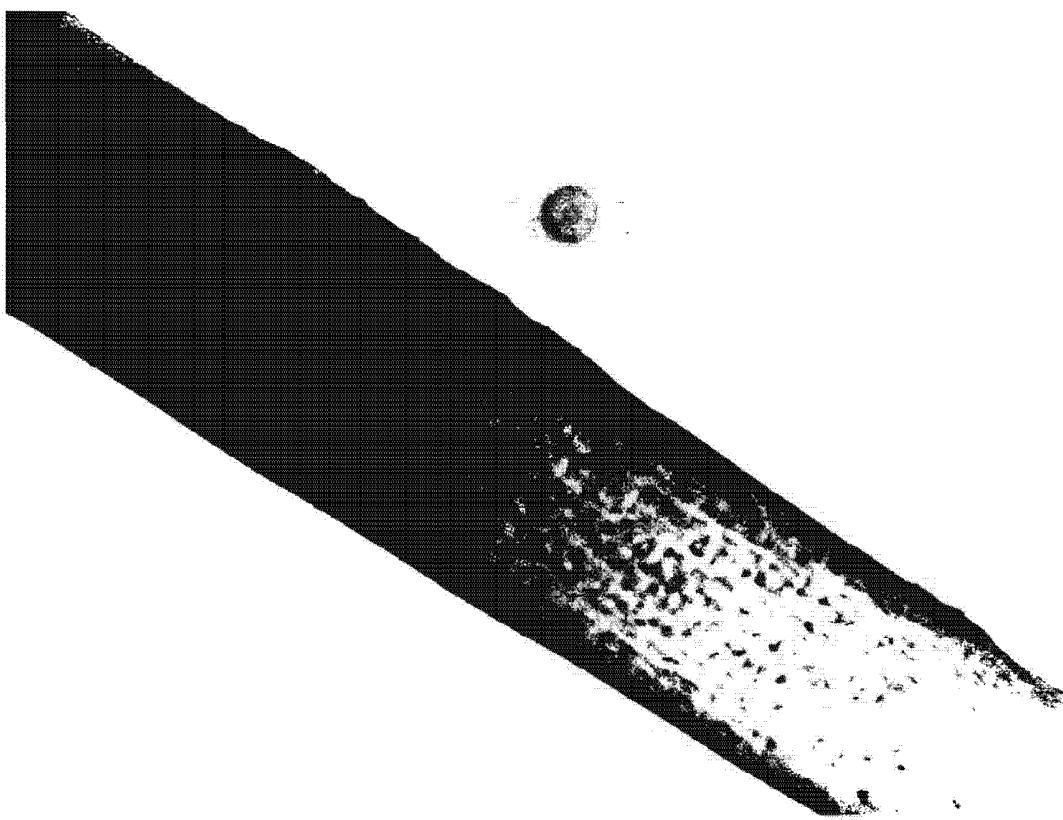

dog 2007 site 3

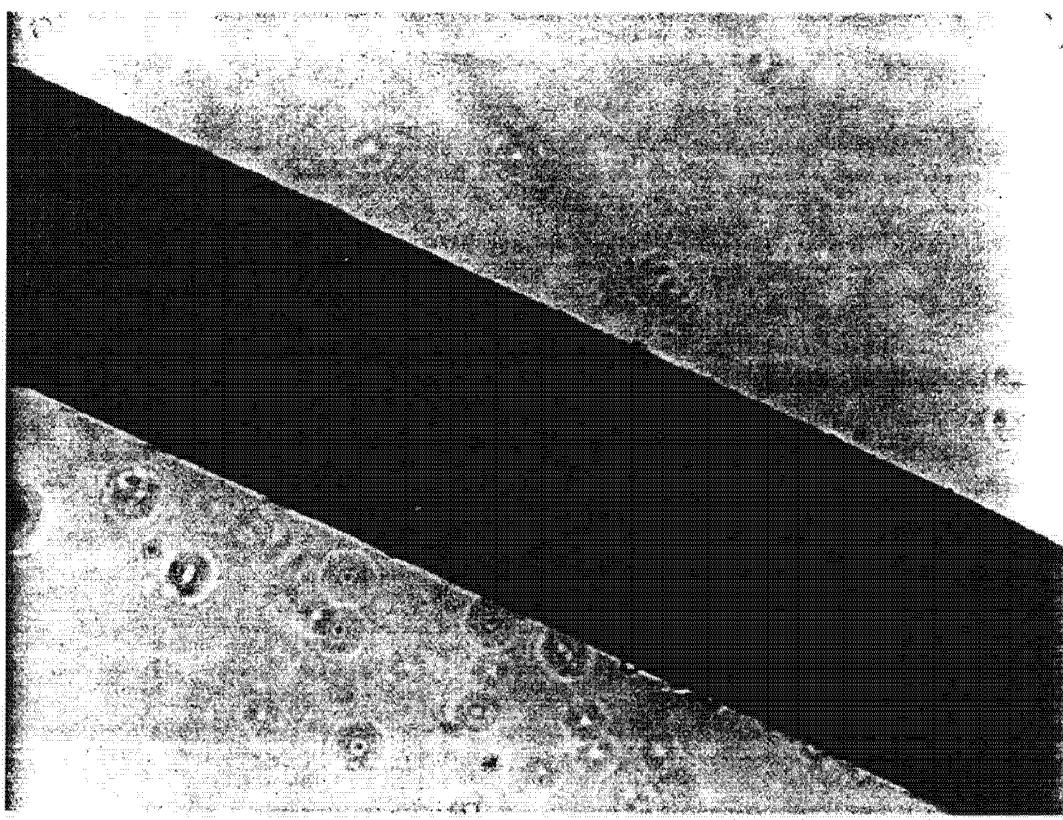



dog 2001 site 2

E.Lindemann 7/15/92 SR# 34097

dog 2008 site 1

Figure 3 7 Year Ethilon Explants

E. Lindemann 7/9/92 SR# 33985*

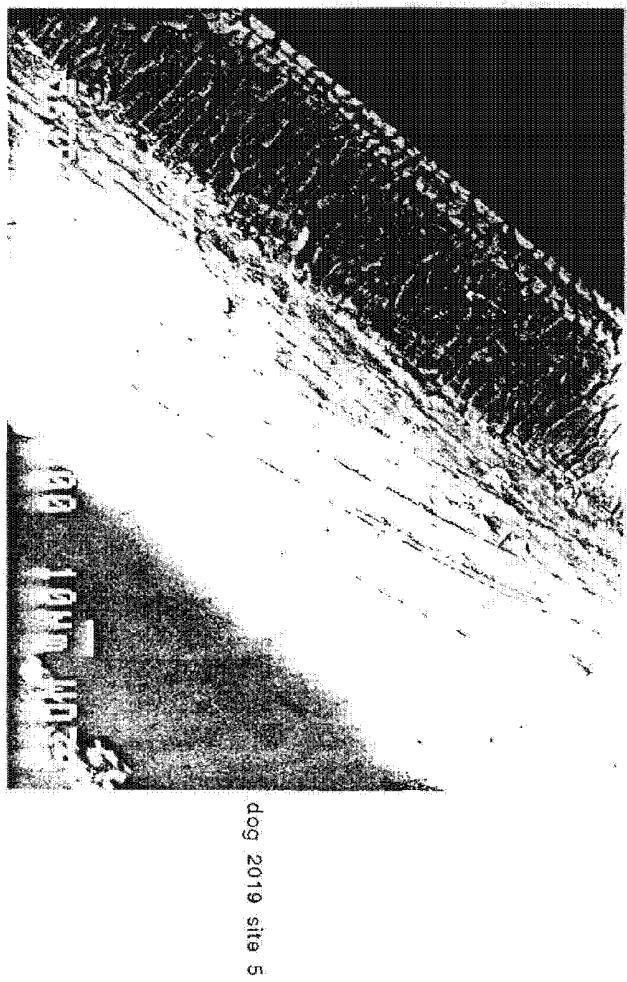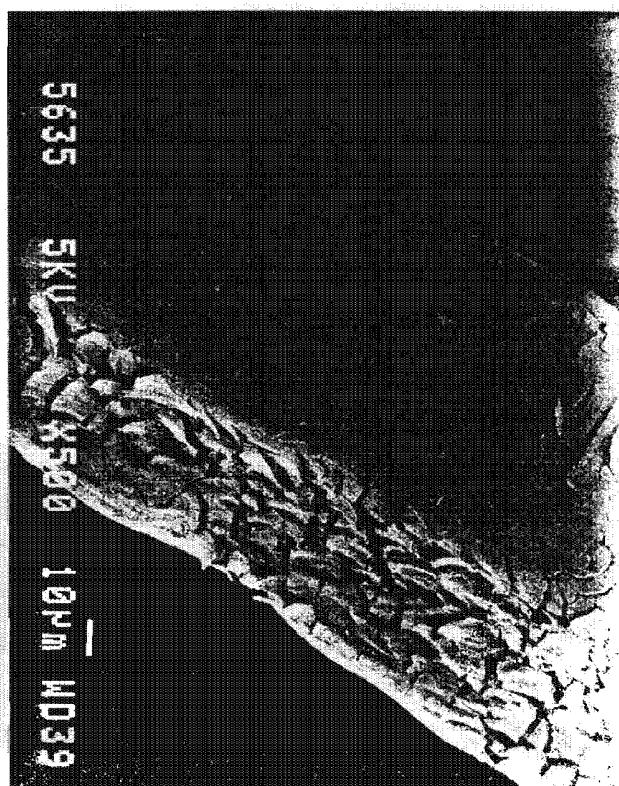

dog 2019 site 5

E. Lindemann 7/9/92 SR# 33985

Figure 4 7 Year Ethilon Explants


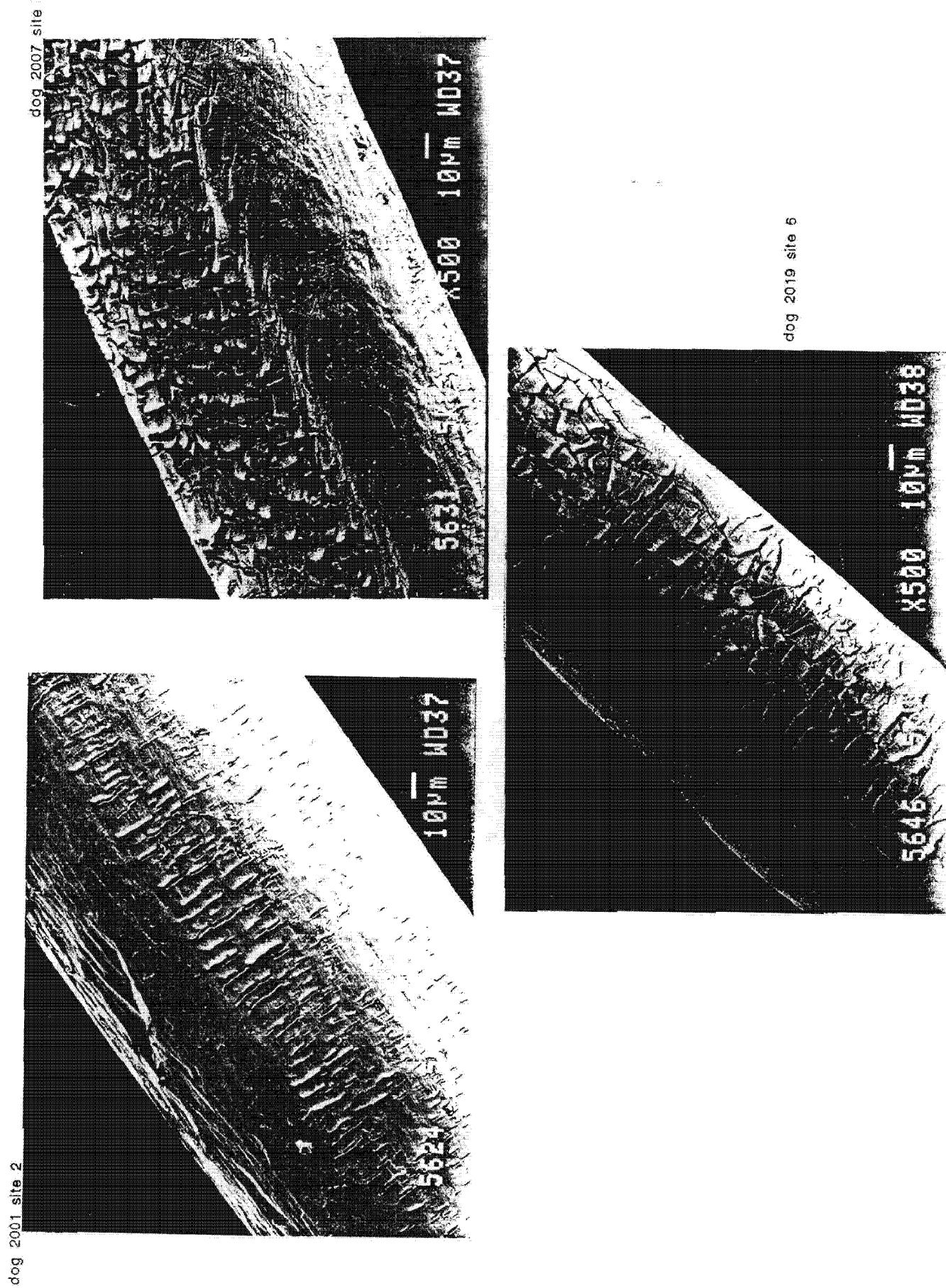


Figure 5A 7 Year PVDF Explants

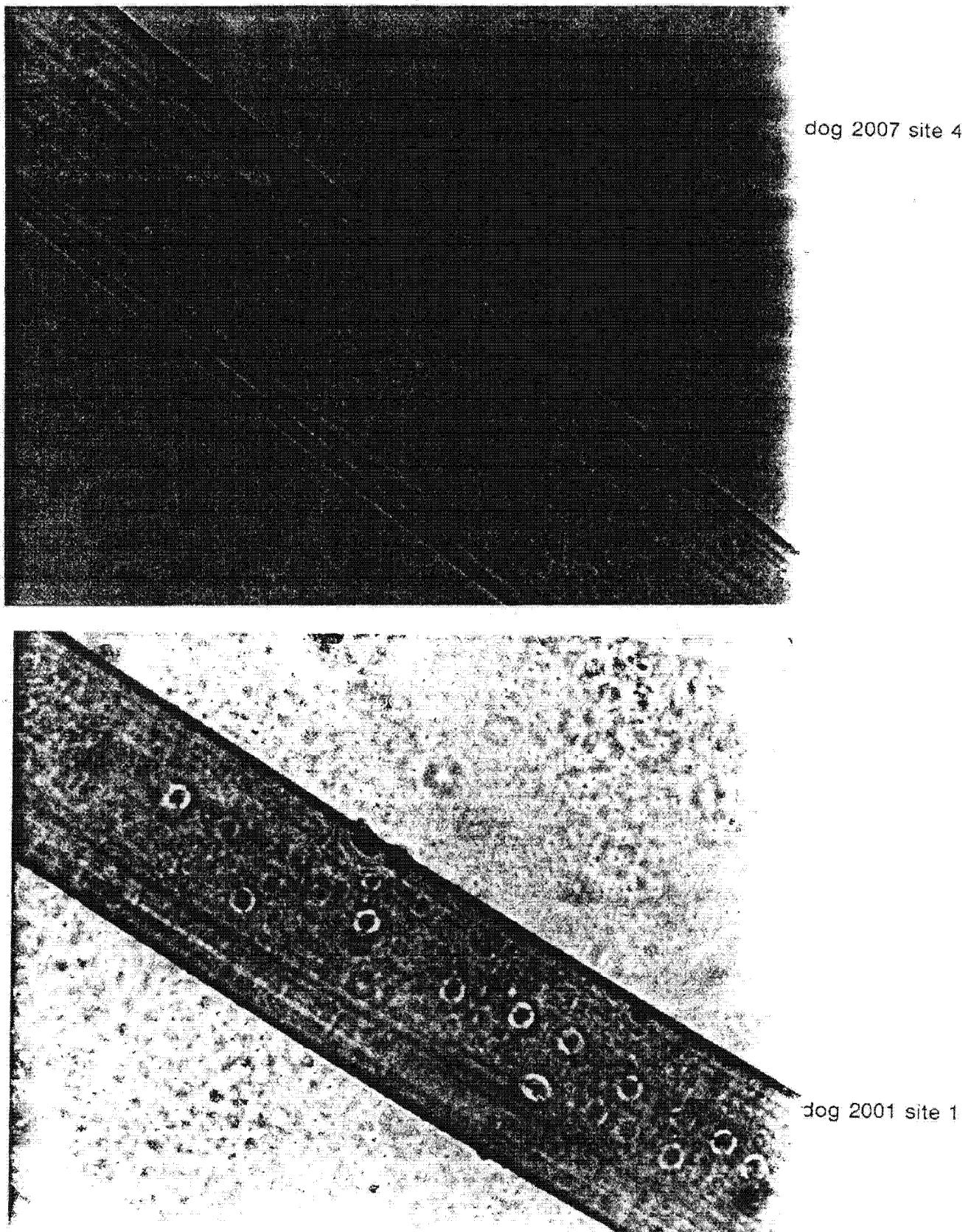


E. Lindemann 7/15/92 SR# 33985


CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690773

Figure 5B 7 Year PVDF Explants

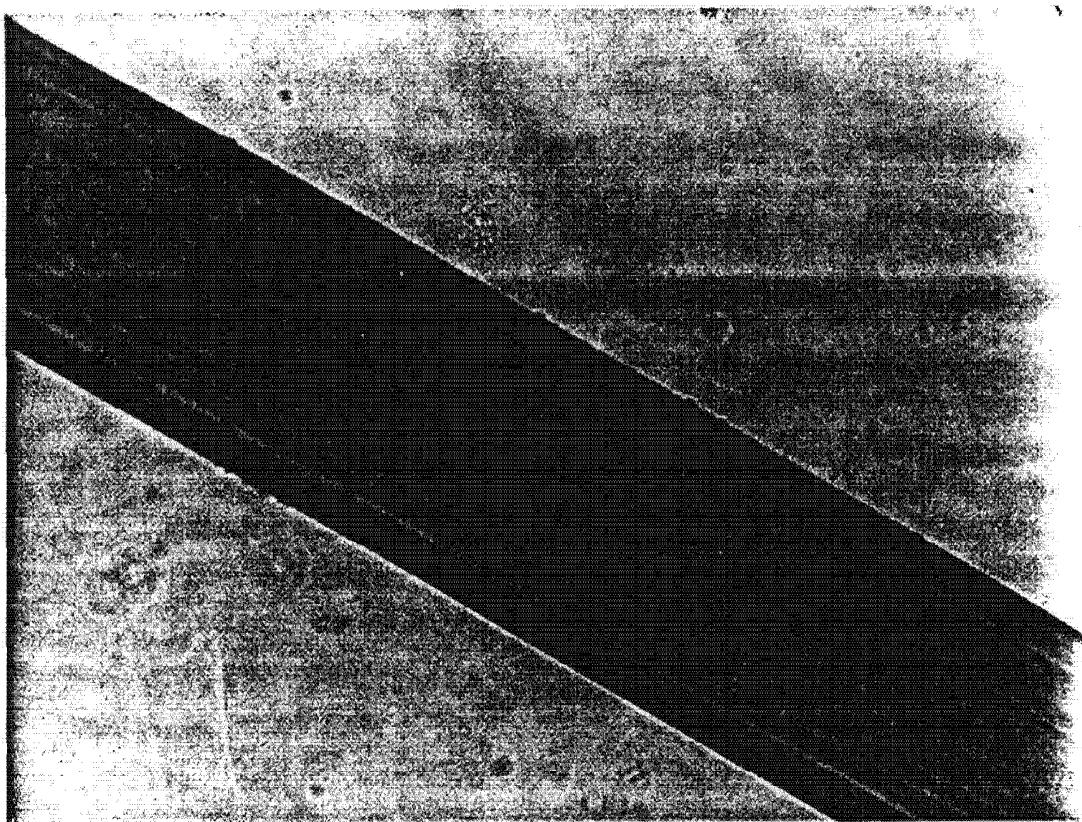

dog 2008 site 4

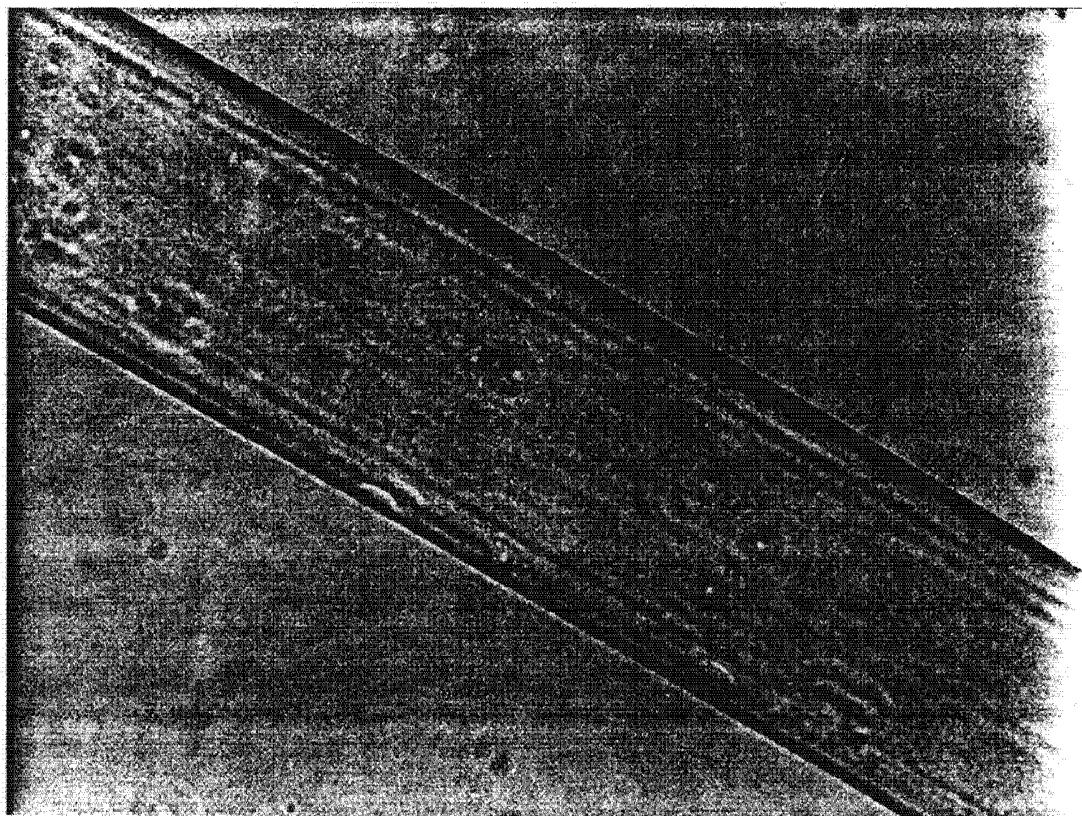



dog 2007 site 5

E. Lindemann 7/15/92 SR# 33985

CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690774

Figure 5C 7 Year PVDF Explants

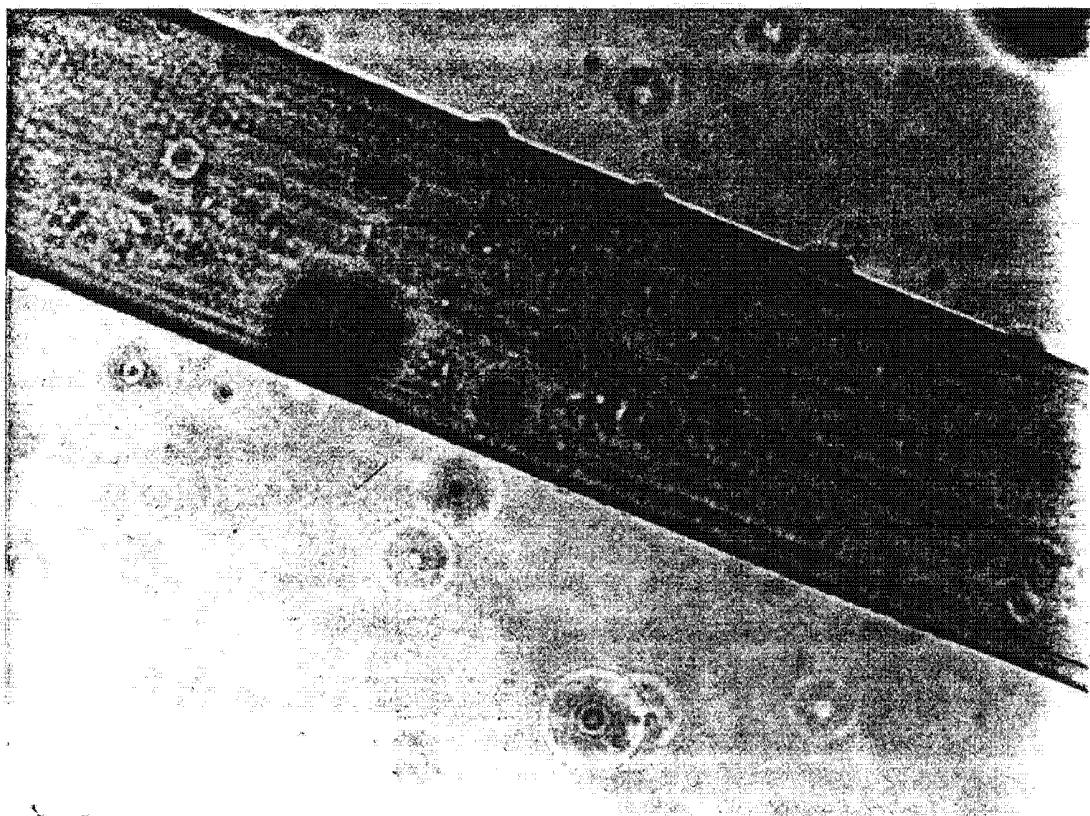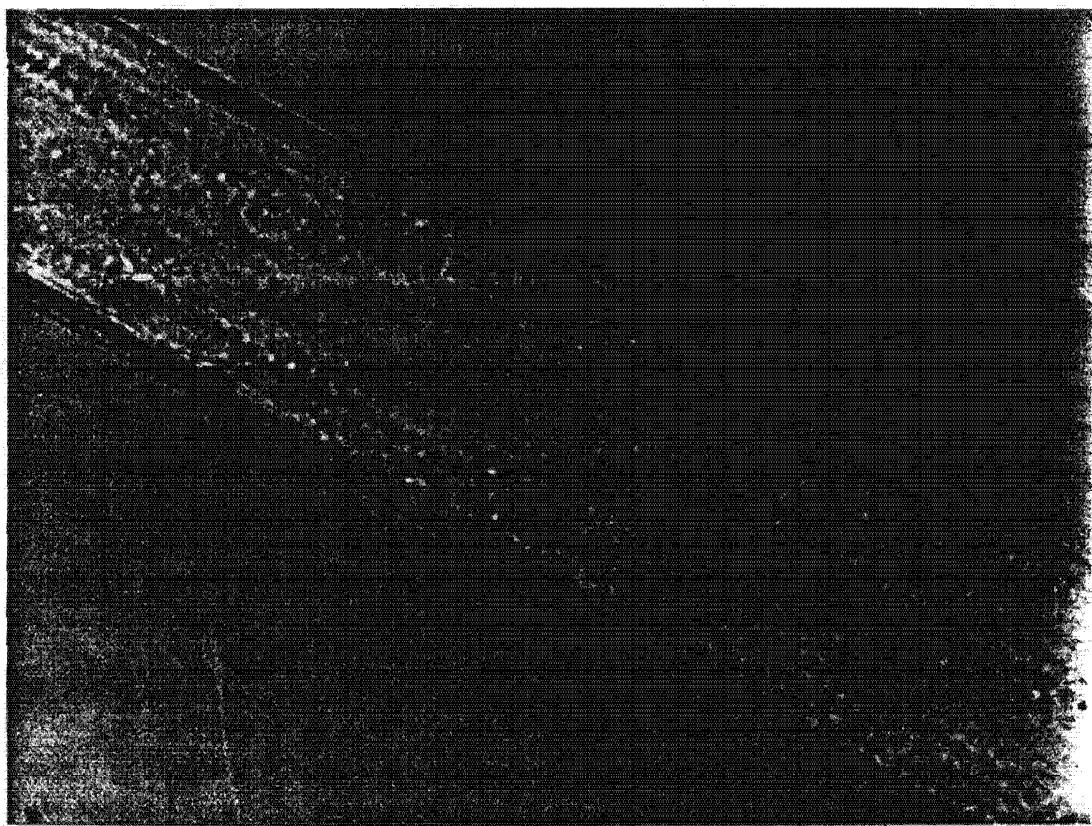


E. Lindemann 7/15/92 SR# 33985

CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690775

E. Lindemann 7/9/92 SR# 33985*

Figure 5 7 Year PVDF Explants


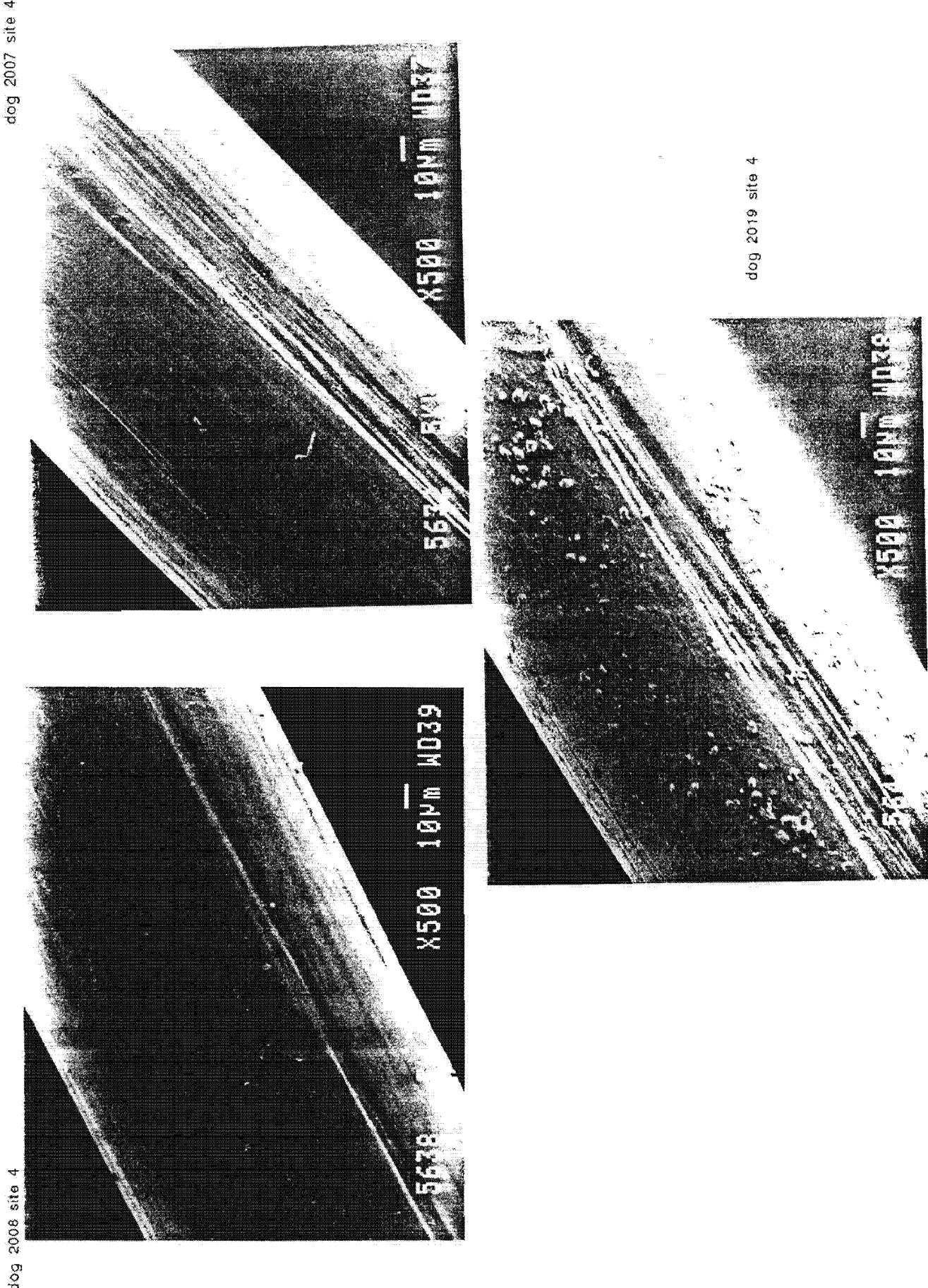


CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690776

E. Lindemann 7/9/92 SR# 33985


Figure 6 7 Year PVDF Explants

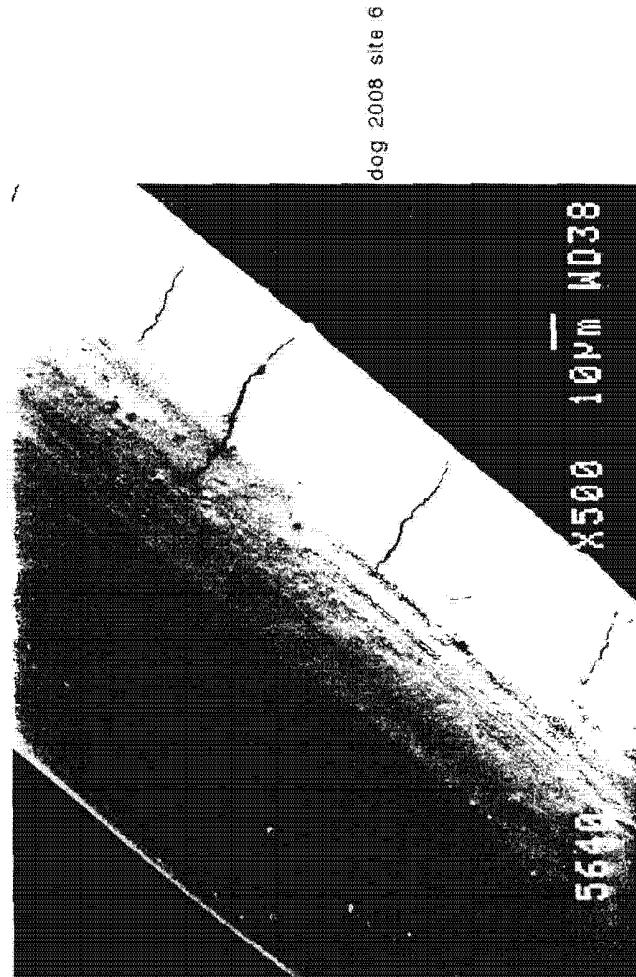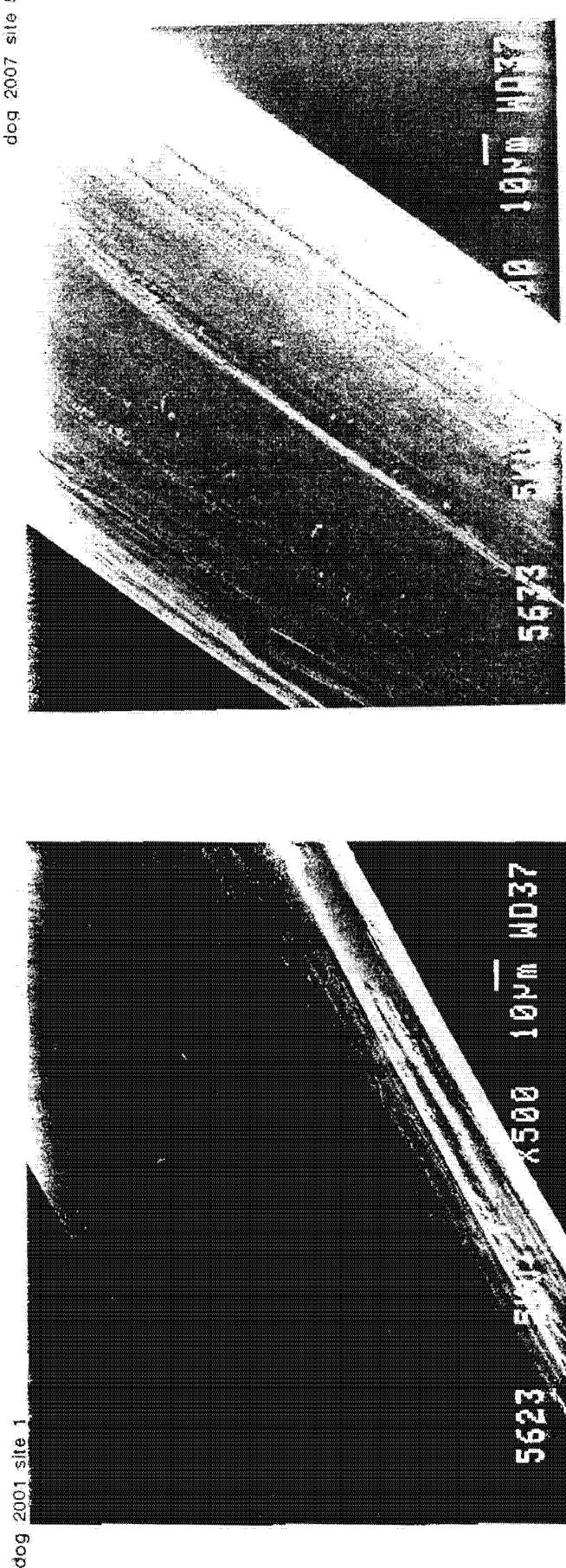


dog 2001 site 1

ETH.MESH.07690777

Figure 7A 7 Year Novafil Explants

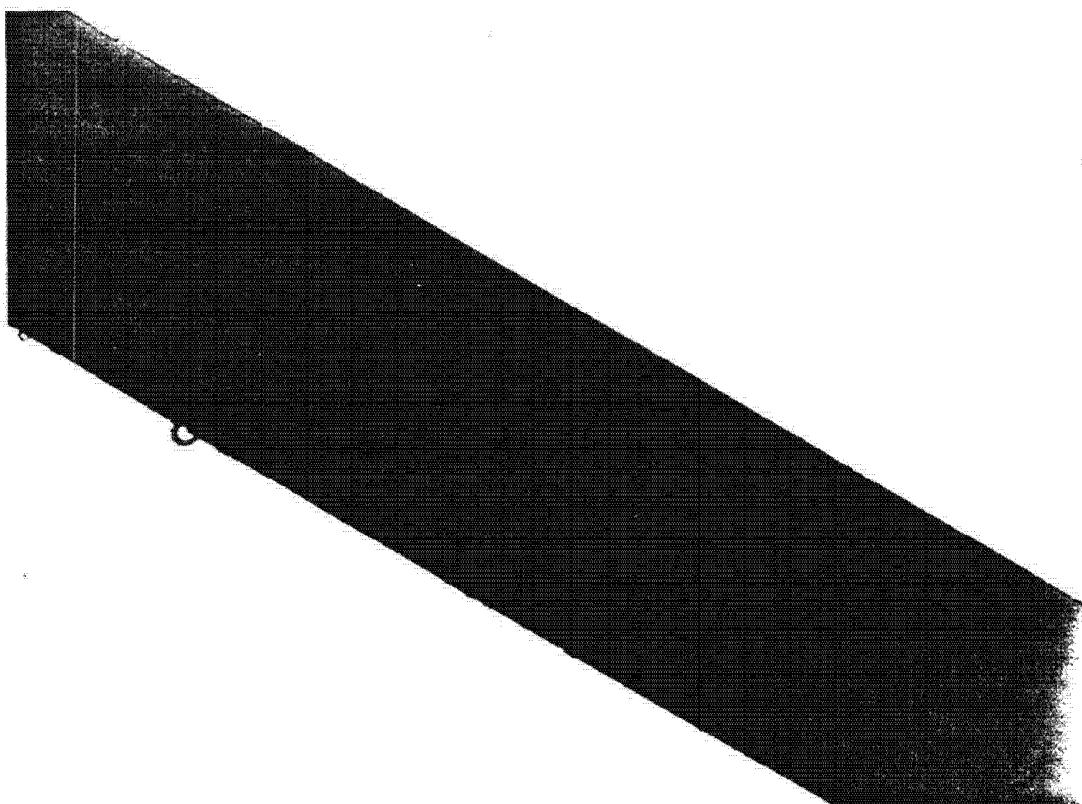

dog 2001 site 6

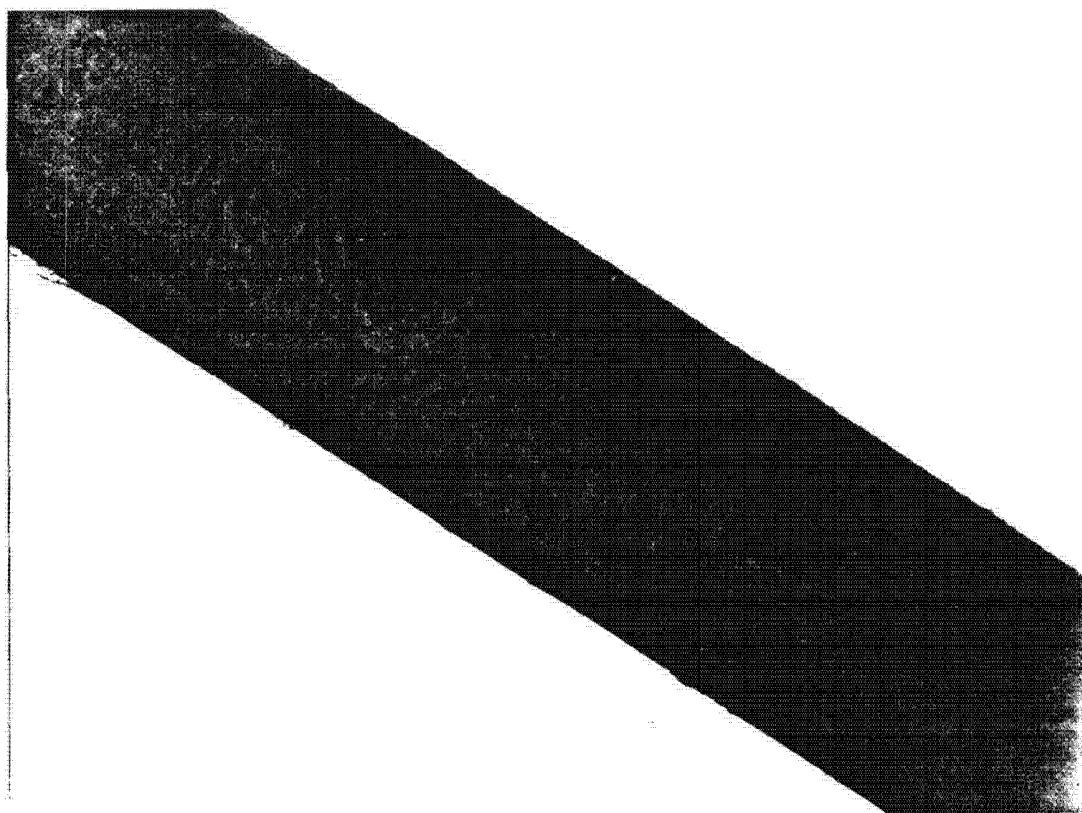


dog 2001 site 3

E. Lindemann 7/15/92 SR# 33985


CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690778

Figure 7B 7 Year Novafil Explants

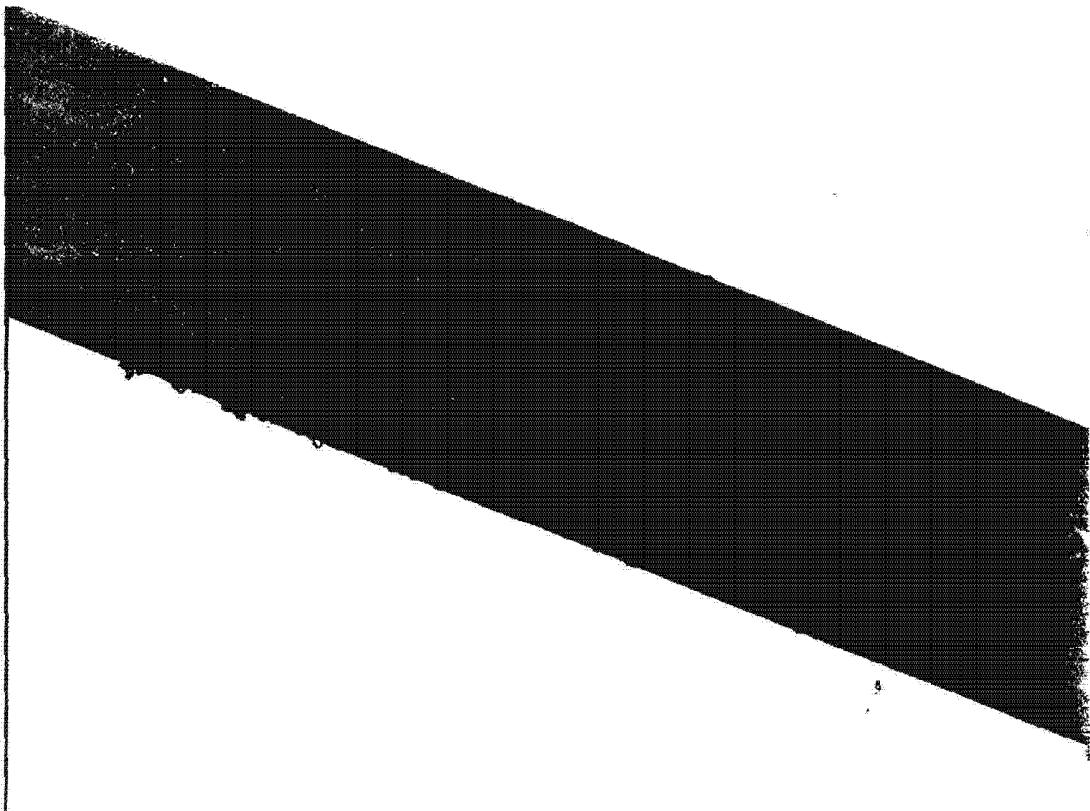

dog 2007 site 2

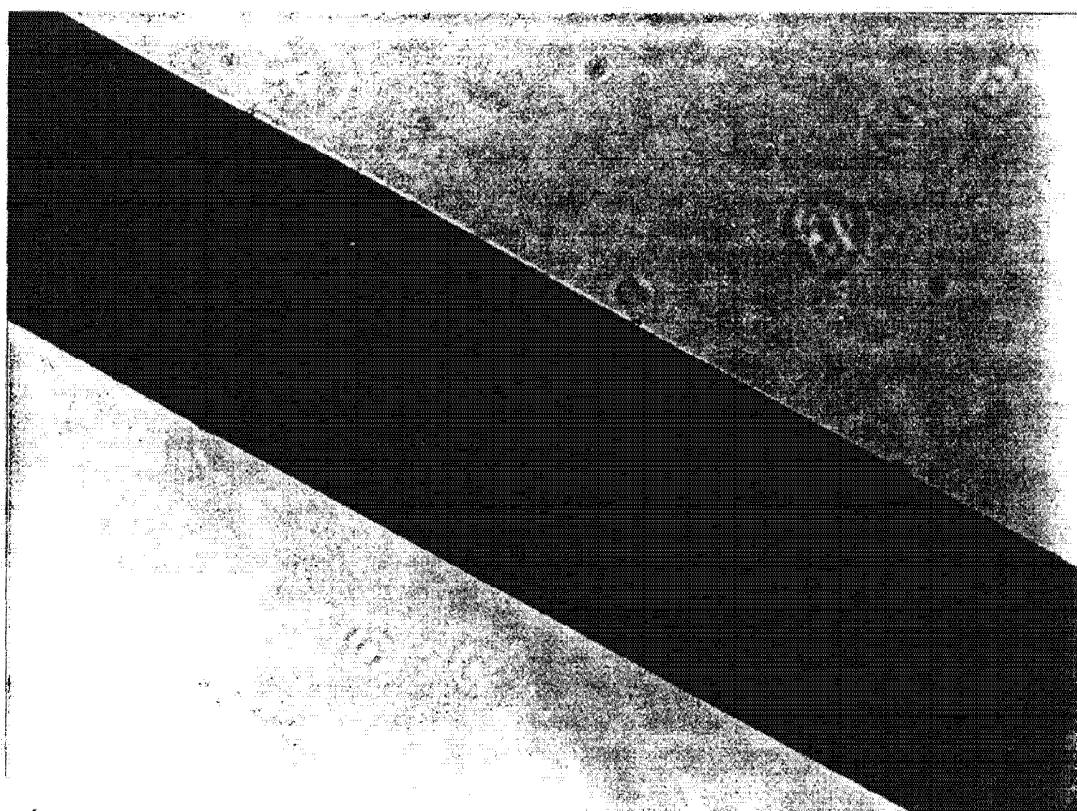

dog 2019 site 1


E. Lindemann 7/15/92 SR# 33985

CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690779

Figure 7C 7 Year Novafil Explants

dog 2008 site 3


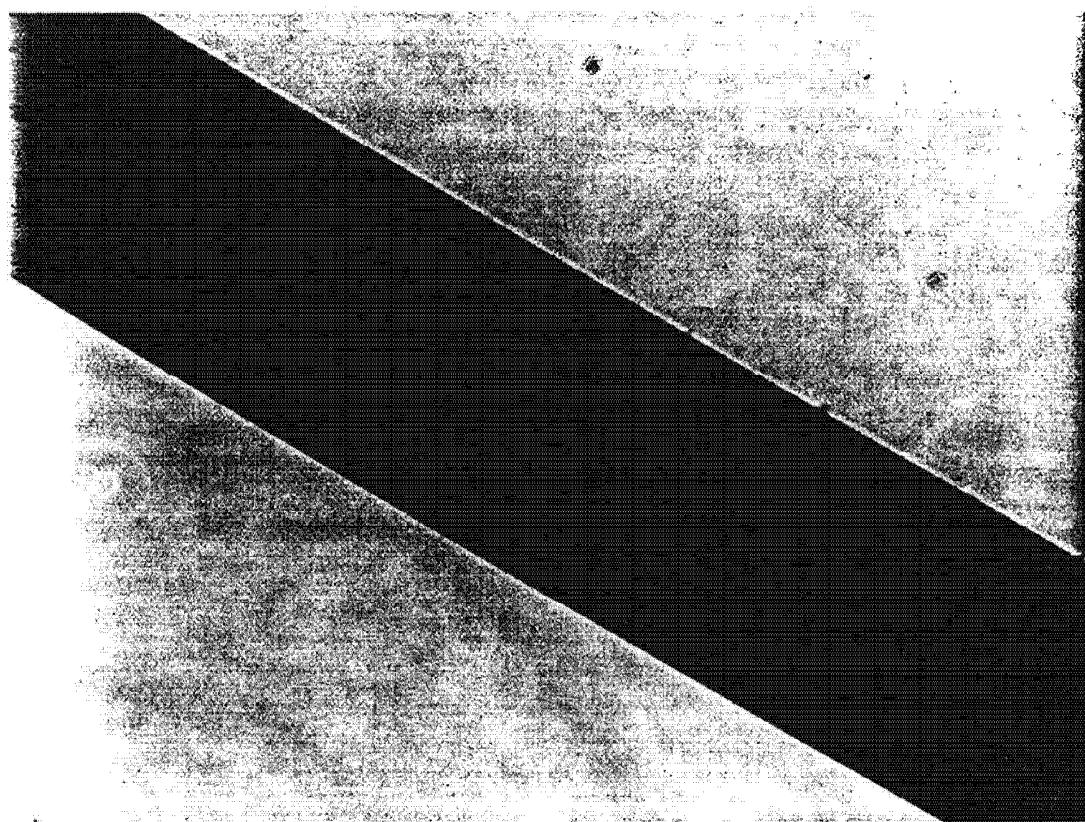


Lindemann 7/15/92 SR# 33985

CONFIDENTIAL
SUBJECT TO STIPULATION AND ORDER OF CONFIDENTIALITY

ETH.MESH.07690780

Figure 7 7 Year Novafil Explants

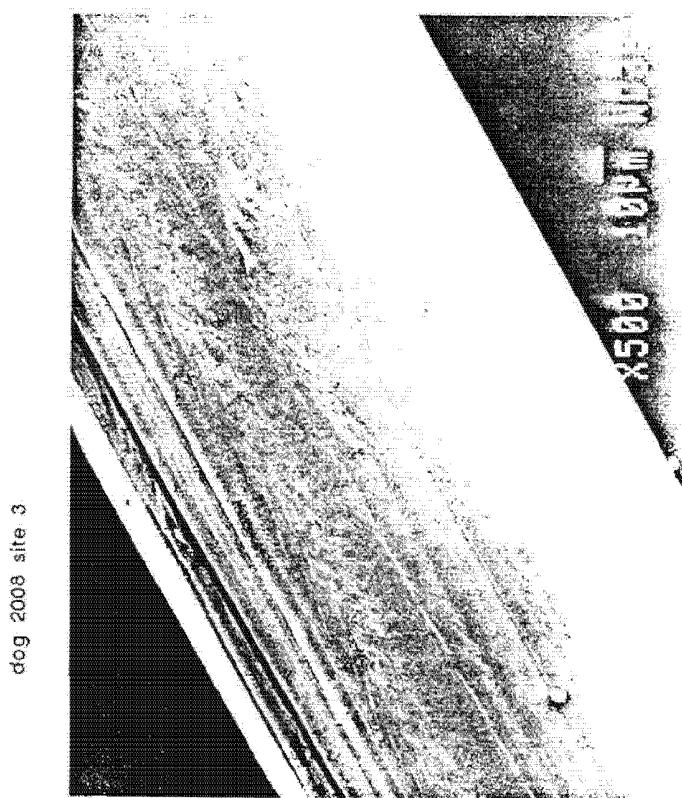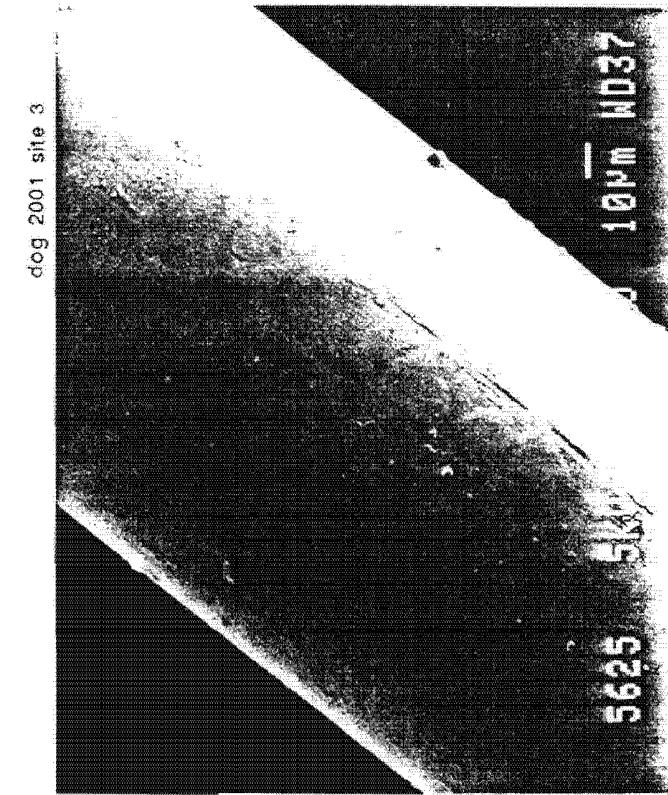

E. Lindemann 7/9/92 SR# 33985

dog 2007 site 2



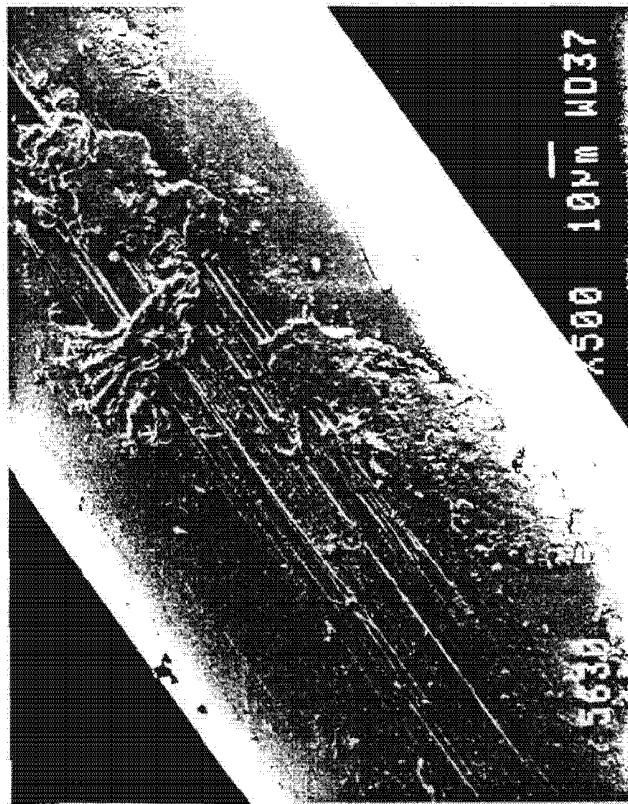


Figure 8 7 Year Novafil Explants

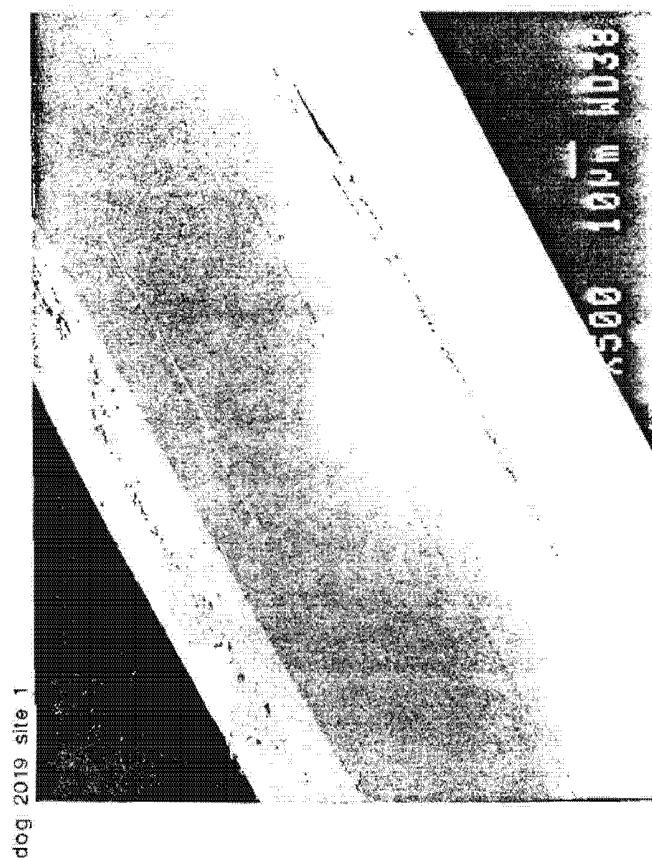

dog 2019 site 1

dog 2001 site 6

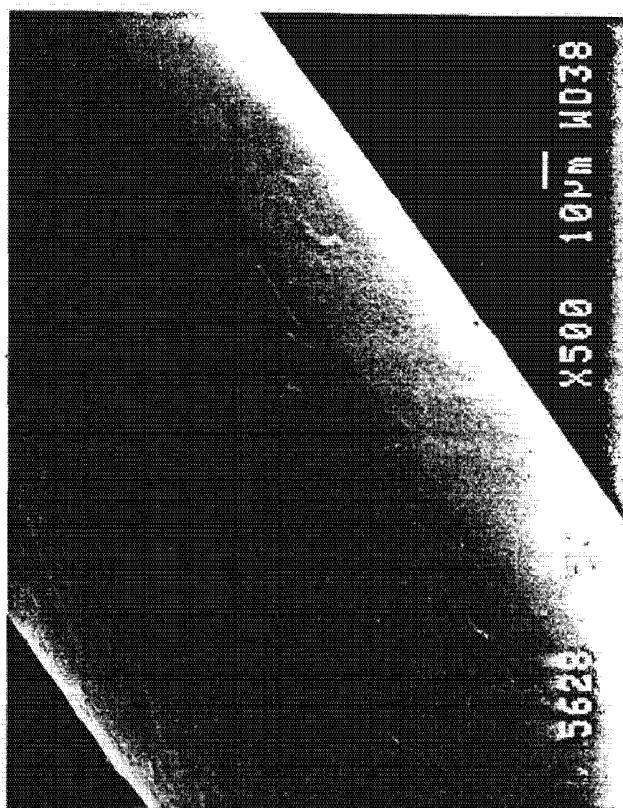

ETHICON, INC.

ANALYTICAL CHEMISTRY DEPARTMENT

SAMPLE CONTROL #

J 3028

SERVICE REQUEST #

3400292

34003

REQUESTOR	DEPARTMENT	EXT.	DATE SUBMITTED	PROJECT NO.	REQUESTOR'S MANAGER AUTHORIZATION
V. AGARWAL	64523	X2205	June 23 '92	16102	B. Matlaga
SIZE 5-0, DOG #2007 7 YEAR EXPLANT					
SAMPLE IDENTIFICATION, SPECIAL STORAGE CONDITIONS, PRECAUTIONS 1. NOVAFIL SITE 2 2. ETHILON SITE 3 3. PVDF (Undyed) SITE 4 4. PVDF (Undyed) SITE 5 5. PROLENE SITE 6 6. PROLENE SITE 1 } 7 Year Explants for [10 year Prolene BSR Study]					
SAMPLE DISPOSITION <input checked="" type="checkbox"/> ATTACHED SENT TO SUPERVISOR/ANALYST O. CONTACT REO. O. RETURN					
ANALYSIS REQUESTED/PURPOSE: IB GC GPC DSC OM % PPM NMR TLC IV TGA HSM O ₂ H ₂ O ETO MS LC XRD TMA SEM DYE MONO UV WET DEN DMA EDXA CRYST COMP OTHER: FOR OM+SEM DATA SEE SR 34070					
NO. SAMPLES <input type="text"/>					

TEST REPORT Samples of each explant were prepared as hot-pressed films by HSM and examined by FT-IR.

SITE 1 - PROLENE : IR spectrum obtained corresponds to polypropylene
 SITE 2 - NOVAFIL : " " "
 SITE 3 - ETHILON : " " "
 SITE 4 - PVDF : " " "
 SITE 5 - Prolene : " " "
 SITE 6 - Prolene : " " "
 SITE 5 - PVDF : " " "

RB 7/8/92

polypropylene
 Novafil (P-1400, 1,4-Butene Tempether
 -polybutoxy ether)
 Nylon 6 (ETHILON).
 polyvinylidene fluoride
 polypropylene
 polypropylene
 polyvinylidene fluoride

Jan 17/86 7/8/92

SAMPLE	IV/d/g	Robin R. Raylana NB 2519-50		
NOVAFIL SITE 2	0.84			
ETHILON SITE 3	1.07			
Insufficient sample for prolene IV				
SAMPLE	IV/d/g	MW	MN	MW/MN
NOVAFIL SITE 2	32,000	18,000	1.8	
ETHILON SITE 3	57,000	27,000	2.1	
Robin R. Raylana NB 2519-52				
Prolene Site 1	322,000	69,000	Current 4/0	324,000 60,000
Prolene Site 6	323,000	63,000		
				6-23-92 PZ

CONCLUSIONS/COMMENTS:
 Comparison of 7 year explants to current 4/0 Prolene sutures indicates no significant degradation

REF. IR file 34003

2562-94

ANALYST Eugene Muse	DATE 10/4/92	ANALYST Robin R. Raylana	DATE 9/21/92	SECTION MANAGER	DATE
ANALYST/SUPERVISOR	DATE	ANALYST/SUPERVISOR Eugene Muse	DATE 9/21/92	DEPT MANAGER	DATE

F140-206 (Rev. 4/85)

ETHICON, INC.

J-2095

SERVICE REQUEST #

33853

ANALYTICAL CHEMISTRY DEPARTMENT

Samples to Gene Mise → Dan Bursley

Pg 1 of 2

REQUESTOR	DEPARTMENT	EXT	DATE SUBMITTED	PROJECT NO.	REQUESTOR'S MANAGER AUTHORIZATION
VISH AGARWAL	64523	X2205	19 th May '92	16102	B. Muttagee

SAMPLE IDENTIFICATION, SPECIAL STORAGE CONDITIONS, PRECAUTIONS						ANALYSIS REQUESTED/PURPOSE:						
<p>Samples from Dog # A95 on ERF 85-219 long-term Prolene study 6 yrs + 10 months</p> <p>NOVAFIL (site 4) NOVAFIL (site 5) Dog Died Prematurely</p> <p>ETHILON (Site 6) PROLENE (Site 3)</p> <p>PVDF (Site 2) ETHILON (Site 1)</p>						IR	UV	GPC	DSC	OM	% H ₂ O	PPM
NMR	TLC	IV	TGA	HSM	0,		ETO					
MS	LC	XRD	TMA	SEM	DYE	MONO	COMP					
U/V	WET	DEN	DMA	EDXA	CRYST							
OTHER:						See attached protocol						
						See SR 33788 for SEM				NO. SAMPLES 6		
SAMPLE DISPOSITION												
<input checked="" type="checkbox"/> ATTACHED	SENT TO	SUPERVISOR/ANALYST										
0 CONTACT REO.	Burkley	McBain / Burkley										
0 RETURN	Muse											

OM + IR Microspectroscopy (for SEM, see § 33788) - IR microspectroscopy was performed on representative samples that exhibited surface cracking. The following samples were examined by IR microspectroscopy:

NOVAFIL (SITE 4) } For these samples, the explant was examined and
ETHILON (SITE 1) } pictured by OM prior to examining by IR microspectroscopy
PROLENE (SITE 3) }

↓ This sample was examined by IR microscopy without any sample preparation - after which it was gold coated and examined by SEM (see SR 33788).

IR spectra obtained from IR microspectroscopy show no differences when comparing cracked ETHION with a "control" region of the same explant surface (a region that is not cracked or discolored). The same can be said for Novafil as well. HOWEVER, expected oxidation absorbances associated with surface cracking would be masked by the ester carbonyl absorbances inherent in these materials. Evidence of oxidation was observed for the PROGEN explant.

All samples examined and commented on as described above were originally submitted with SR 33788.

Janet B.L. 5/22/92

5-1992 R

15

CONCLUSIONS/COMMENTS:	REF.				
Samples must be refrigerated until tested. Thank you. <i>Approved</i>					
ANALYST <i>Dawn J. B.</i>	DATE <i>5/22/92</i>	ANALYST	DATE	SECTION MANAGER	DATE
ANALYST/SUPERVISOR	DATE	ANALYST/SUPERVISOR	DATE	DEPT. MANAGER <i>Ed M. Lewis</i>	DATE <i>9/1/92</i>

SAMPLE CONTROL #
J2095

ETHICON, INC.

ANALYTICAL CHEMISTRY DEPARTMENT

PAGE 2 OF 2

SERVICE REQUEST NO.	ANALYSIS REQUIRED	ANALYTICAL SUPERVISOR
33853	GPC, IR Identity	

TEST REPORT

The explanted samples were prepared as hot pressed films for IR identity. From the IR spectra generated, the explant samples are identified as:

{

ERF 85-214
Dog 1995

SITE 1 : ETHILON (Nylon 6)
 SITE 2 : PVDF (polyvinylidene fluoride)
 SITE 3 : Prolene (polypropylene)
 SITE 4 : Novafil (polybutylethylene terephthalate-polybutoxypolyether)
 SITE 5 : Novafil
 SITE 6 : ETHILON (Nylon 6)

Janet F. B. 6/15/92
(R file 33853)

SAMPLE	IV/dlg
Novofil (Site 4)	0.73
Novofil (Site 5)	0.82
ETHILON (Site 1)	1.25
ETHILON (Site 6)	1.24

Robin R. Ray/lambda NB 2519-31

Insufficient sample for proline IV's

SAMPLE	MW	MN	MW/MN
Novofil (site 4)	32,000	18,000	1.8
Novofil (site 5)	33,600	18,500	2.0 ^{over} 1.8
ETHILON (site 1)	62,000	30,000	2.1
ETHILON (site 6)	61,000	31,000	2.0

Robin R. Ray/lambda NB 2519-34

Although there was insufficient sample to run IV
GPC of Dog # 1995 Site 3 was compared to a ^{current} 4% Proline suture

	Mw	MN	N
Dog #1995 Site 3	327,000	59,000	
Current 4% Prolene	324,000	60,000	

Results indicate no degradation has taken place EPM 2562-94.

ANALYST SIGNATURE	DATE	REFERENCE	
Robin R. Ray/lambda	9/21/92		

J6542

ANALYTICAL CHEMISTRY DEPARTMENT

34180

REQUESTOR V. AGARWAL	DEPARTMENT 64523	EXT. X2205	DATE SUBMITTED July 27 '92	PROJECT NO. 16102	REQUESTOR'S MANAGER AUTHORIZATION
(DOG 2019) (7 YEAR EXPLANT)					
SAMPLE IDENTIFICATION, SPECIAL STORAGE CONDITIONS, PRECAUTIONS			ANALYSIS REQUESTED/PURPOSE		
1. Prolene Site 2 2. ETHILON Site 5 3. ETHILON Site 6 4. NOVAFIL Site 1 5. PVDF Site 4			Size 5-0 samples DOG # 2019 7 year explant Samples for 10 year Prolene BSR study		
			IR NMR MS U/V	GC TLC LC WET	GPC TV XRD DEN
			DSC TGA TMA DMA	OM HGM SEM EDXA	% O, H ₂ O DYE CRYST
			OTHER	SEM + OM DATA Documented ON SR 34097	
SAMPLE DISPOSITION <input checked="" type="checkbox"/> ATTACHED <input type="checkbox"/> CONTACT REQ. <input type="checkbox"/> RETURN	SENT TO JFB	SUPERVISOR/ANALYST MUSE BURKLEY LITERATOR	NO SAMPLES		

TEST REPORT Samples prepared as hot-pressed films using the hot stage:

Dog 2019, SITE 1 : IR spectrum corresponds to Novafil.

SITE 2 : IR spectrum corresponds to PROLENE.

SITE 3 : " " " PROLENE

SITE 4 : " " " PVDF

SITE 5 : " " " ETHILON

SITE 6 : " " " ETHILON

insufficient sample for prolene IV

Janet F B 8/5/92

SAMPLE IV/dig

ETHILON SITES 1.00

ETHILON SITE 6 1.16

NOVAFIL SITE 1 0.87

Robin Rhey/anal NB 2519-69

Insufficient sample for prolene IV

SAMPLE MW MN MW/MN

ETHILON SITES 72,000 21,000 3.4

ETHILON SITE 6 74,000 22,000 3.3

NOVAFIL SITE 1 32,000 13,000 2.4

Robin Rhey/anal NB 2519-71

Current Prolene 4/6 324,000 60,000

Dog # 2019 Site 3 331,000 64,000

Dog # 2019 Site 2 332,000 57,000

7-28-92

CONCLUSIONS/COMMENTS

Comparison of 7 year explants to current prolene indicate no molecular weight degradation

REF IR file 34180
2562-94

ANALYST Janet F B	DATE 8/5/92	ANALYST Robin R. Rhey/anal	DATE 9/21/92	SECTION MANAGER	DATE
ANALYST/SUPERVISOR Eugene Muse	DATE 10/4/92	ANALYST/SUPERVISOR Eugene Muse	DATE 9/21/92	DEPT. MANAGER	DATE

EHA 206 (Rev 4/85)

ANALYTICAL CHEMISTRY DEPT COPY

SAMPLE CONTROL #
I 2261

SERVICE REQUEST #
34066

ETHICON, INC.
ANALYTICAL CHEMISTRY DEPARTMENT

REQUESTOR	DEPARTMENT	EXT.	DATE SUBMITTED	PROJECT NO.	REQUESTOR'S MANAGER AUTHORIZATION
V. AGARWAL	64523	2205	July 2'92	16102	

SAMPLE IDENTIFICATION, SPECIAL STORAGE CONDITIONS, PRECAUTIONS 7 Year Explant Samples from 10-Year Protocol Size 5-0 DOG # 2008 BSR Study						ANALYSIS REQUESTED/PURPOSE	
① Prolene Site #2 ② PVDF Site #6 ③ ETHILON Site #5 ④ PVDF Site #4 ⑤ NOVAFIL Site #3 ⑥ ETHILON Site #1						IR GC GPC DSC OM % : PPM NMR TLC IV TGA HSM O ₂ H ₂ O ETO MS LC XRD TMA SEM DYE MONO UIV WET DEN DMA EDXA CRYST COMP OTHER:	
SAMPLE DISPOSITION		SENT TO	SUPERVISOR/ANALYST				NO. SAMPLES
<input checked="" type="checkbox"/> ATTACHED <input type="checkbox"/> CONTACT REQ <input type="checkbox"/> RETURN			<i>McWayne Baptist Memphis</i>				

1. S. Please keep the samples refrigerated until tested

TEST REPORT FOR IR Identity, a sample of each explant was hot-pressed into a film on the hot stage, followed by IR examination.

SITE 1 - ETHILON: IR spectrum obtained corresponds to Nylon 6

SITE 2 - PROLENE: " " " " polypropylene

SITE 3 - Novafil: " " " " polybutylene terephthalate-polybutylene

SITE 4 - PVDF: " " " " polyvinylidene fluoride

SITE 5 - ETHILON: " " " " Nylon 6

SITE 6 - PVDF: " " " " polyvinylidene fluoride

all IR spectra verify the identity of the sutures implanted.

7/9/92 *Janet F. B.*
MW

SAMPLE IV/d19

ETHILON #1 # 0.95

Dog# 2008 Site 2 322,000

53,000

ETHILON #5 1.28

Current Proleve #0 324,000

60,000

NOVAFIL #3 0.79

Robin R Ray/anh NO 2519-50

2562-94

Insufficient sample for prolene IV

SAMPLE MW MN MW/MN

ETHILON site #1 53,000 27,000 2.0

ETHILON site #5 59,000 28,000 2.1

Novafil site #3 32,000 18,000 1.8

Robin R Ray/anh NB 2519-52

7-2-92

CONCLUSIONS/COMMENTS:

Comparison of current prolene 4/0 suture indicates no significant degradation of 7 yr explant

REF IR fl 34066

ANALYST <i>JFB</i>	DATE 7/11/92	ANALYST <i>Robin R. Ray/anh</i>	DATE 9/21/92	SECTION MANAGER	DATE
ANALYST/SUPERVISOR <i>Eugene Muse</i>	DATE 10/9/92	ANALYST/SUPERVISOR <i>Eugene Muse</i>	DATE 9/21/92	DEPT. MANAGER	DATE

E140-206 Rev 4/85

ANALYTICAL CHEMISTRY DEPT. COPY

ACCESSION 85-219
PROJECT NO. 16102

EXPLANTATION PROCEDURES and SAMPLE DISTRIBUTION

Notify the following people of upcoming explant dates:

Nancy Myirski, x2743: Microscopic inspection - someone from her group will come to inspect the sutures under the dissecting scope under the hood. Samples should be placed after dissection from dog into saline-moistened paper towels labelled with the ERF acc. no., dog no., site no., suture type and date.
(Ann Leibold was inspector @ 2 yr. time period, 6/87)

Frank Schiller, x3040: SEM - An Analytical Chemistry Service Request form must be filled out and accompany each set of samples.
Put the sample control number on the top left corner of the sample label. Make one copy for our file and one to send with the samples. Mail original to Dr. A. Melvegar. Label samples the same as for above.

Implantation (Stef or Dan?)

Kevin Sullivan, x2997: Instron - Submit samples after the microscopic inspection, while moist. Fragments are saved in their respective towels for next tests. Refrigerate if there will be a delay between inspection and instroning.

Gene Muse, x3046: Molecular weight - Deliver moist suture fragments after Instroning. After testing he will deliver samples to Dan.

Dan Burkley, x3048: I.R. - Receives samples from Gene. Will discard samples when testing is completed.

Explant samples in consecutive order. Dissect both LC100's (dorsal and ventral) from surrounding connective tissue, carefully stripping tissue from the suture surface. Cut one of the LC100's off the sutures at the clip and gently pull the suture bundle through the tissue by gripping the remaining LC100. When free of tissue, moisten with saline and separate one strand from the bundle. Place this strand into a large (15 ml) red-top tube filled with sterile water and labelled as described above. The other 5 strands per bundle are placed in moistened paper towels labelled as described above. The single sutures in tubes are submitted for SEM and the remaining strands are inspected microscopically and tested on the Instron, etc. as described above.

* 2 request forms per dog: - 1 for SEM + OM + IR microscopy for the samples going to J. Melvegar
Have JOP sign each & send to A. Melvegar.
- 1 for IR, GPC + IV for the fragments going to G. Muse + D. Burkley